



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

TSY  
3.10.2010

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002759 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3/8/2010 Cal 1 0.100 AC Cal 2 0.100 AC

Between the dates of - and - the instrument was operating properly and accurately.

The instrument was taken out of service for repair on 10/22/2009.  
The instrument was returned to service on 3/8/2010.

Sarah Cione

Sarah L. Cione


3/10/10

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 10th day of March, 20 10.

Susan A. Hitchye

Notary Public

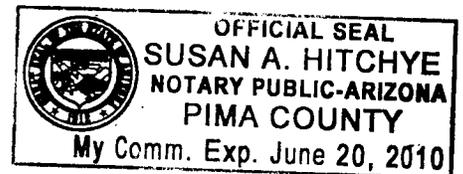

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

1 of 16  
SC

INTOXILYZER SERVICE RECORD

QAS: Sarah Cione

BH  
3-10-2010

Intoxilyzer Serial Number: 80-002759

Date Returned to Service: 3/8/2010

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD QUALITY ASSURANCE PROCEDURE

TSY  
3-10-2010

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 3/8/2010 TIME 1031

INTOXILYZER SERIAL # 80-002759 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.

(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.

(✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed

(✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.

(✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 100 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Sarah Cione

3 of 16 sc

INTOXILYZER 8000  
Location:  
Serial Number : 80-002759  
03/08/2010

10:49:07

Standard Lot#: 826001I  
Last Changed By: S. CIONE #50616

DN

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:49:53
Diagnostic Test	Pass	10:50:21
Air Blank	0.000	10:50:48
0.100 Cal Check	0.101	10:51:07
Air Blank	0.000	10:51:36
Subject Test	0.000	10:52:14
Air Blank	0.000	10:52:41
Five Minute Wait	ABT*	
Air Blank	0.000	10:53:29

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

4 of 16 SC

INTOXILYZER 8000  
Location:  
Serial Number : 80-002759  
03/08/2010

10:57:22

AS

---

Standard Lot#: 826001I  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	RFI*	10:58:12
Air Blank	IPS**	10:58:22

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

5 of 16 sc

INTOXILYZER 8000

Location:

Serial Number : 80-002759

03/08/2010

10:58:50

AS4

Standard Lot#: 826001I

Last Changed By: S. CIONE

#50616

QAS: S. CIONE

#50616

TUCSON PD

Operator: S. CIONE

#50616

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:59:50
Diagnostic Test	Pass	11:00:18
Air Blank	0.000	11:00:45
0.100 Cal Check	0.100	11:01:04
Air Blank	0.000	11:01:33
Subject Test	0.102*	11:04:41
Air Blank	0.000	11:05:13
Air Blank	0.000	11:05:41
Subject Test	XXX**	11:06:14
Air Blank	0.000	11:06:48

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000

Location:

Serial Number: 80-002759

03/08/2010

11:07:50

-----  
Standard Lot#: 826001I

Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	11:08:18
Cal Check	0.100	11:08:36
Air Blank	0.000	11:09:06

  
\_\_\_\_\_  
Operator's Signature

6 of 16 sc

BSH

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

7 of 16 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

TS4  
3.10.2010

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 3/8/2010 TIME 1031  
Intoxilyzer Serial # 80-002759 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location:

Serial Number: 80-002759

03/08/2010

10:31:46

-----  
Standard Lot#: 826001I

Last Changed By: S. CIONE #50616

-----  
QAS: S. CIONE #50616

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:32:13
Cal Check	0.100	10:32:32
Air Blank	0.000	10:33:01

*Sarah Cione*

Operator's Signature

8 of 16 sc

124



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

9 of 16 SC

INTOXILYZER SERVICE RECORD  
FOR  
CALIBRATION GAS TANK CHANGE

TSY  
3-10-2010

Date: 3/8/2010

QAS: Sarah Cione

Intoxilyzer Serial Number: 80-002759

New Lot Number: 826001 I

New Expiration Date: 9/24/2010

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."  
Select "D" for dry gas and enter the following information:  
Enter Target Value: 0.100  
Enter Control Lot #: From the tank label  
Enter Exp Date: From the tank label (mm/dd/yyyy)  
Edit info Y/N? Y or N  
Update Std Info? Y ENTER  
"Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "Control Testing."  
Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0.100 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

10 of 16 sc

INTOXILYZER 8000

Location:

Serial Number: 80-002759

03/08/2010

10:30:03

-----  
Standard Lot#: 826001I

Last Changed By: S. CIONE #50616

-----  
QAS: S. CIONE #50616

TUCSON PD  
-----

TS4

Test	g/210L	Time
Air Blank	0.000	10:30:31
Cal Check	0.100	10:30:50
Air Blank	0.000	10:31:19

Sarah Cione  
Operator's Signature



12 of 16 sc

154

# Certificate of Calibration

## **Intoxilyzer**® 8000

*This is to certify that the calibration of **Intoxilyzer**® serial number 80-002759, manufactured by CMI Inc., a subsidiary of MPD, Inc. of Owensboro, Kentucky, was tested and found to be in compliance with the National Highway Traffic Safety Administration Standard for Devices to Measure Breath Alcohol (F.R., No. 179 48705-48710 Sept. 17, 1993) and with International Recommendation OIML R126 Edition 1998 (E). Calibration solutions are traceable to NIST standard material 1828.*

Date 12/11/09

Signed

*Kristi Robinson*

Technician

**CMI** INC.

316 East 9<sup>th</sup> Street  
Owensboro, KY 42303  
USA

Part No. 650519 CMI 7/9/09

Pursuant to Section 28.1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

13 of 16 SC

INTOXILYZER SERVICE RECORD

TS4  
3-10-2010

Date: 12/23/09

QAS: Sarah Cione

Intoxilyzer Serial Number: 80-002759

Reason For Service: Received from CMI on 12/23/09, but instrument still has Diagnostic failures (Voltage/Current)

Date Out of Service: -

Calibration Check completed prior to removal from service

Date: -

SQAP completed prior to removal from service

Date: -

Corrective Action:

Sent to CMI, Inc.

sent back to CMI



# Return Material Authorization Form

14 of 16 SC

154  
3.10.2010

This form MUST be completed and enclosed with item(s) being returned. Failure to complete and return this form may cause delays in crediting / repairing.

1. Item(s) Returned:  Instrument  Supplies  Other  
(Note: please ship items in their original shipping container(s) or a similar protective box.)

2. Instrument Model: Intox 8000 Serial Number 80-002759

3. Bill To Address: Finance / Accounts Payable Ship To Address: (NO PO BOXES)  
P.O. Box 27450 Tucson Police Department  
TUCSON, AZ 85726-7450 Crime Lab  
270 S. Stone Ave  
TUCSON, AZ 85701-1917

Customer # 105075

Name Terry Gallegos

Phone: (520) 837-7677 Fax: (520) 791-4642

Email: Terry.Gallegos@tucsonaz.gov

4. Reason for Return: Note - If you are sending an item for repair, please give a detailed description of the problem. Please list any special instructions that you may have concerning this return. Note - Returns for credit may be subject to a 15% restocking charge.

Received from CMI 12/23/09, but instrument having same problem (fails  
Diagnostics - voltage/current)

5. Purchase Order Number (attach a copy of P.O. if applicable) 20669

6. I Authorize All Repairs:  Yes  No

7. I Authorize Repairs Up To: \$ \_\_\_\_\_

8. I Require An Estimate Regardless of Cost  Yes  No

Note: An evaluation fee will apply to estimates that are not repaired.

9. Authorized By:

Terry Gallegos  
Name (Please Print)

[Signature]  
Signature

12.28.2009  
Date

Pursuant to Section 28.1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.



16 of 16 SC

TS4  
3.10.2010

# Certificate of Calibration

## **Intoxilyzer**® 8000

*This is to certify that the calibration of **Intoxilyzer**® serial number 80-002759, manufactured by CMI Inc., a subsidiary of MPD, Inc. of Owensboro, Kentucky, was tested and found to be in compliance with the National Highway Traffic Safety Administration Standard for Devices to Measure Breath Alcohol (F.R., No. 179 48705-48710 Sept. 17, 1993) and with International Recommendation OIML R126 Edition 1998 (E). Calibration solutions are traceable to NIST standard material 1828.*

Date 2/24/10 Signed Kristi Robinson  
Technician

**CMI** INC.

316 East 9<sup>th</sup> Street  
Owensboro, KY 42303  
USA

Part No. 650519 CMI 7/9/09

Pursuant to Section 28.1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.