



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

TSY
5-10-2010

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002759 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 05/06/10 Cal 1 0.102 AC Cal 2 0.099 AC

Between the dates of 04/08/10 and 05/06/10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Cheryl K. Straub

Cheryl K. Straub


05/06/10

ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 6th, day of May, 2010.

A.M. Hernandez

Notary Public

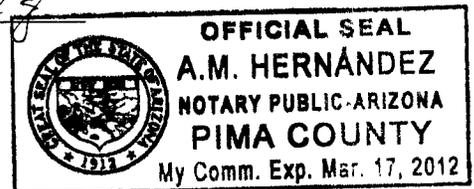

*Includes documents specified in A.A.C. R13-10-104.A.5

187
CS

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

TSY
5-10-2010

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Cheryl STRAUB AGENCY Tucson Police Department

DATE 05/06/10 TIME 10:34

INTOXILYZER SERIAL # 80-002759 LOCATION CRIME LAB

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Cheryl Straub

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

277
05

TS4

INTOXILYZER 8000
Location:
Serial Number : 80-002759
05/06/2010 10:36:50

Standard Lot#: 826001I
Last Changed By: S. CIONE #50616

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:37:41
Diagnostic Test	Pass	10:38:08
Air Blank	0.000	10:38:35
0.100 Cal Check	0.101	10:38:54
Air Blank	0.000	10:39:24
Subject Test	0.000	10:39:59
Air Blank	ABT*	10:40:06
Air Blank	0.000	10:40:35

*Sequence Aborted

Not a Successfully Completed Test Sequence

377
CS

124

INTOXILYZER 8000

Location:

Serial Number : 80-002759

05/06/2010

10:46:30

Standard Lot#: 826001I

Last Changed By: S. CIONE

#50616

QAS: C. STRAUB

#50772

TUCSON PD

Operator: C. STRAUB

#50772

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	10:47:37
Air Blank	IPS**	10:47:45

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

487
CS

BS4

INTOXILYZER 8000

Location:

Serial Number : 80-002759

05/06/2010

10:48:10

Standard Lot#: 826001I

Last Changed By: S. CIONE

#50616

QAS: C. STRAUB

#50772

TUCSON PD

Operator: C. STRAUB

#50772

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:49:02
Diagnostic Test	Pass	10:49:29
Air Blank	0.000	10:49:56
0.100 Cal Check	0.101	10:50:15
Air Blank	0.000	10:50:45
Subject Test	0.590*	10:53:53
Air Blank	0.000	10:54:28
Air Blank	0.000	10:54:56
Subject Test	XXX**	10:55:19
Air Blank	0.000	10:55:52

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

577
08 48708
08

INTOXILYZER 8000

Location:

Serial Number: 80-002759

05/06/2010

10:56:14

Standard Lot#: 826001I

Last Changed By: S. CIONE #50616

QAS: C. STRAUB #50772

TUCSON PD

Test g/210L Time

Air Blank 0.000 10:56:42
Cal Check 0.099 10:57:00
Air Blank 0.000 10:57:30



Operator's Signature

637
CS

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

T54
5-10-2010

QA Specialist Cheryl Straub AGENCY Tucson Police Department
DATE 05/06/10 TIME 10:34
Intoxilyzer Serial # 80-002759 Location CRIME LAB

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution N/A AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 102 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Cheryl Straub

DPS Form Exh G-2 (Rev. 05-01)

