



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

TS4
6.8.2010

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002759 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6/3/10 Cal 1 0.100 AC Cal 2 0.100 AC

Between the dates of 5/6/10 and 6/3/10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


6/3/10

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 3rd, day of June, 2010.

A.M. Hernandez

Notary Public

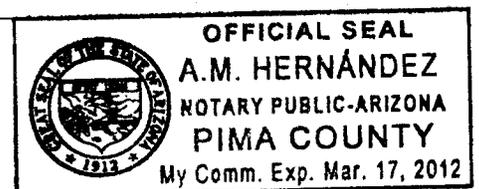

*Includes documents specified in A.A.C. R13-10-104.A.5

1 of 7 sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TS4
6-8-2010

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 6/3/10 TIME 1015

INTOXILYZER SERIAL # 80-002759 LOCATION Team 4

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

INTOXILYZER 8000
Location:
Serial Number : 80-002759
06/03/2010

10:16:35

2 of 7 sc

TS4

Standard Lot#: 719919
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|------------------|--------|----------|
| Air Blank | 0.000 | 10:17:21 |
| Diagnostic Test | Pass | 10:17:48 |
| Air Blank | 0.000 | 10:18:15 |
| 0.100 Cal Check | 0.101 | 10:18:34 |
| Air Blank | 0.000 | 10:19:04 |
| Subject Test | 0.000 | 10:19:22 |
| Air Blank | 0.000 | 10:19:50 |
| Five Minute Wait | ABT* | |
| Air Blank | 0.000 | 10:20:24 |

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location:

Serial Number : 80-002759

06/03/2010

10:20:30

3 of 7 sc

TS4

Standard Lot#: 719919

Last Changed By: C. CONLEY

#34517

QAS: S. CIONE

#50616

TUCSON PD

Operator: S. CIONE

#50616

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | RFI* | 10:20:58 |
| Air Blank | IPS** | 10:21:05 |

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location:
Serial Number : 80-002759
06/03/2010 10:21:11

4 of 7 sc
TS4

Standard Lot#: 719919
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 10:21:53 |
| Diagnostic Test | Pass | 10:22:20 |
| Air Blank | 0.000 | 10:22:48 |
| 0.100 Cal Check | 0.101 | 10:23:06 |
| Air Blank | 0.000 | 10:23:36 |
| Subject Test | 0.089* | 10:26:44 |
| Air Blank | 0.000 | 10:27:16 |
| Air Blank | 0.000 | 10:27:43 |
| Subject Test | XXX** | 10:27:57 |
| Air Blank | 0.000 | 10:28:30 |

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location:

Serial Number: 80-002759

06/03/2010

10:29:23

5 of 7 SL

TSY

Standard Lot#: 719919

Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616

TUCSON PD

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | 0.000 | 10:29:51 |
| Cal Check | 0.100 | 10:30:09 |
| Air Blank | 0.000 | 10:30:39 |

Sarah Cione

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6 of 7 sc
TS4
6-8-2010

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department

DATE 6/3/10 TIME 1015

Intoxilyzer Serial # 80-002759 Location Team 4

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location:

Serial Number: 80-002759

06/03/2010

10:15:13

7 of 7 sc
TSh

Standard Lot#: 719919

Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616

TUCSON PD

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | 0.000 | 10:15:41 |
| Cal Check | 0.100 | 10:15:59 |
| Air Blank | 0.000 | 10:16:29 |

Sarah Cione

Operator's Signature



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

1 of 2
3C

INTOXILYZER SERVICE RECORD
FOR
CALIBRATION GAS TANK CHANGE

TS4
6.8.2010

Date: 5-7-10
QAS: Clare Conley
Intoxilyzer Serial Number: 80-002759
New Lot Number: 719919
New Expiration Date: 02/02/2012

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."
Select "D" for dry gas and enter the following information:
 - Enter Target Value: 0.100
 - Enter Control Lot #: From the tank label
 - Enter Exp Date: From the tank label (mm/dd/yyyy)
 - Edit info Y/N? Y or N
 - Update Std Info? Y ENTER
 - "Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."
Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 101 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove and sign printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

2 of 2 sc

BY

INTOXILYZER 8000

Location:

Serial Number: 80-002759

05/07/2010

08:08:14

Standard Lot#: 719919

Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517

TUCSON PD

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | 0.000 | 08:08:42 |
| Cal Check | 0.101 | 08:09:01 |
| Air Blank | 0.000 | 08:09:30 |

Clare Conley

Operator's Signature