



CERTIFICATION by CUSTODIAN of RECORDS  
 Tucson Police Department Crime Laboratory  
 270 S. Stone Avenue Tucson, AZ 85701-1917

*By  
2-10-2011*

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002759 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 02/10/2011 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of 01/13/2011 and 02/10/2011 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.  
 The instrument was returned to service on N/A.

*Cheryl K. Straub*

Cheryl K. Straub


02/10/2011

ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 10, day of February, 20 11.

*A.M. Hernandez*

Notary Public

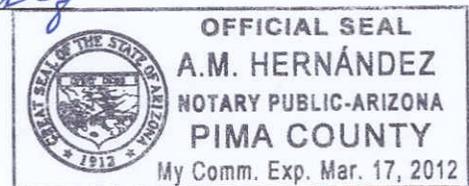

\*Includes documents specified in A.A.C. R13-10-104.A.5

187cs

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

TS4  
2.10.2011

QA SPECIALIST Cheryl STRAUB AGENCY Tucson Police Department

DATE 02/10/2011 TIME 06:00

INTOXILYZER SERIAL # 80-002759 LOCATION TEAM 4

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 099 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Cheryl Straub

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

277 CS

INTOXILYZER 8000

Location:

Serial Number : 80-002759

02/10/2011

06:02:04

BY

Standard Lot#: 719919

Last Changed By: C. CONLEY

#34517

OAS: C. STRAUB

#50772

TUCSON PD

Operator: C. STRAUB

#50772

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	06:03:02
Diagnostic Test	Pass	06:03:30
Air Blank	0.000	06:03:57
0.100 Cal Check	0.100	06:04:16
Air Blank	0.000	06:04:46
Subject Test	0.000	06:05:04
Air Blank	ABT*	06:05:11
Air Blank	0.000	06:05:40

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

377 CS

INTOXILYZER 8000

Location:

Serial Number : 80-002759

02/10/2011

06:05:45

BH

Standard Lot#: 719919

Last Changed By: C. CONLEY

#34517

QAS: C. STRAUB

#50772

TUCSON PD

Operator: C. STRAUB

#50772

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	06:06:17
Air Blank	IPS**	06:06:27

\*RFI Detect

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1527, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

477  
CS

Bh

INTOXILYZER 8000

Location:

Serial Number : 80-002759

02/10/2011

06:06:32

Standard Lot#: 719919

Last Changed By: C. CONLEY

#34517

QAS: C. STRAUB

#50772

TUCSON PD

Operator: C. STRAUB

#50772

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	06:07:22
Diagnostic Test	Pass	06:07:49
Air Blank	0.000	06:08:16
0.100 Cal Check	0.101	06:08:35
Air Blank	0.000	06:09:04
Subject Test	0.071*	06:12:13
Air Blank	0.000	06:12:45
Air Blank	0.000	06:13:12
Subject Test	XXX**	06:13:42
Air Blank	0.000	06:14:16

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

Pursuant to section 26-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

577 CS

154

INTOXILYZER 8000

Location:

Serial Number: 80-002759

02/10/2011

06:14:48

-----  
Standard Lot#: 719919

Last Changed By: C. CONLEY #34517  
-----

QAS: C. STRAUB #50772

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	06:15:15
Cal Check	0.099	06:15:34
Air Blank	0.000	06:16:04



Operator's Signature

EXHIBIT G-2

687cs

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

TS4  
2-10-2011

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Cheryl STRAUB AGENCY Tucson Police Department  
DATE 02/10/2011 TIME 06:00  
Intoxilyzer Serial # 80-002759 Location TEAM 4

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution N/A AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 100 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Cheryl Straub

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

787  
CS

INTOXILYZER 8000

Location:

Serial Number: 80-002759

02/10/2011

06:00:39

154


-----  
Standard Lot#: 719919

Last Changed By: C. CONLEY #34517

-----  
QAS: C. STRAUB #50772

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	06:01:07
Cal Check	0.100	06:01:26
Air Blank	0.000	06:01:55

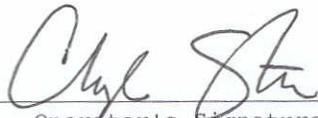  
Operator's Signature