



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002759 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4
3-15-2011

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 03/10/2011 Cal 1 0.107 AC Cal 2 0.099 AC

Between the dates of 02/10/2011 and 03/10/2011 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Cheryl K. Straub

Cheryl K. Straub


03/15/2011

ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 15, day of March, 2011.

A.M. Hernandez

Notary Public

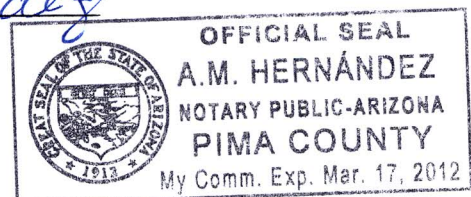

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

187
CSTHIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

TS4
3-15-2011QA SPECIALIST Cheryl Straub AGENCY Tucson Police DepartmentDATE 03/10/2011 TIME 05:39INTOXILYZER SERIAL # 80-002759 LOCATION TEAM 4

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
-
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0.
- 000
- AC.
-
- (✓) 2. Error recognition logic system functioning.
-
- Not a Successfully Completed Test Sequence printed
-
- (✓) 3. Proper sample recognition system.
-
- Not a Successfully Completed Test Sequence printed
-
- Deficient sample printed.
-
- (✓) 4. Standard Calibration Check standard 0.
- 100
- AC
-
- Result: 0.
- 099
- AC

Instrument is operating properly and accurately. YES ✓ NO

COMMENTS

SIGNATURE

Cheryl Straub

277
CS

INTOXILYZER 8000

Location:

Serial Number : 80-002759

03/10/2011

05:40:37

154

Standard Lot#: 719919

Last Changed By: C. CONLEY

#34517

QAS: C. STRAUB

#50772

TUCSON PD

Operator: C. STRAUB

#50772

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	05:41:35
Diagnostic Test	Pass	05:42:02
Air Blank	0.000	05:42:30
0.100 Cal Check	0.099	05:42:48
Air Blank	0.000	05:43:18
Subject Test	0.000	05:43:35
Air Blank	ABT*	05:43:42
Air Blank	0.000	05:44:11

*Sequence Aborted

387
CS

INTOXILYZER 8000

Location:

Serial Number : 80-002759

03/10/2011

05:44:16

154

Standard Lot#: 719919

Last Changed By: C. CONLEY

#34517

QAS: C. STRAUB

#50772

TUCSON PD

Operator: C. STRAUB

#50772

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	05:44:47
Air Blank	IPS**	05:44:55

*RFI Detect

**Improper Sample

477
CS

INTOXILYZER 8000

Location:

Serial Number : 80-002759

03/10/2011

05:45:02

ph

Standard Lot#: 719919

Last Changed By: C. CONLEY

#34517

QAS: C. STRAUB

#50772

TUCSON PD

Operator: C. STRAUB

#50772

TUCSON PD

Subject: SOAP, SOAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	05:45:52
Diagnostic Test	Pass	05:46:20
Air Blank	0.000	05:46:47
0.100 Cal Check	0.100	05:47:06
Air Blank	0.000	05:47:35
Subject Test	0.076*	05:50:44
Air Blank	0.000	05:51:16
Air Blank	0.000	05:51:43
Subject Test	XXX**	05:52:04
Air Blank	0.000	05:52:38

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

587
CS

INTOXILYZER 8000

Location:

Serial Number: 80-002759

03/10/2011

05:53:12

154

Standard Lot#: 719919

Last Changed By: C. CONLEY #34517

QAS: C. STRAUB #50772

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	05:53:39
Cal Check	0.099	05:53:58
Air Blank	0.000	05:54:28



Operator's Signature

EXHIBIT G-2

697
CS

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TS4
3.15.2011

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Cheryl STRAUB AGENCY Tucson Police Department
DATE 03/10/2011 TIME 05:39
Intoxilyzer Serial # 80-002759 Location TEAM 4

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution N/A AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 107 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Cheryl Straub

DPS Form Exh G-2 (Rev. 05-01)

777
cs

INTOXILYZER 8000

Location:

Serial Number: 80-002759

03/10/2011

05:39:11

54


Standard Lot#: 719919

Last Changed By: C. CONLEY #34517

QAS: C. STRAUB #50772

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	05:39:39
Cal Check	0.107	05:39:58
Air Blank	0.000	05:40:27

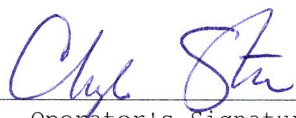

Operator's Signature