



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002759 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4
4.13.2011

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 04/07/11 Cal 1 0.099 AC Cal 2 0.097 AC

Between the dates of 03/10/2011 and 04/07/2011 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Cheryl K. Straub

Cheryl K. Straub


04/07/2011

ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 7th day of April, 2011.

A.M. Hernandez

Notary Public

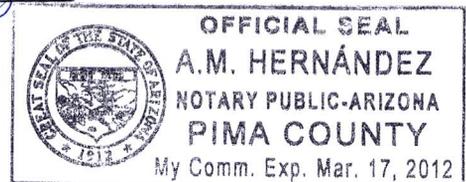

*Includes documents specified in A.A.C. R13-10-104.A.5

197
CS

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD QUALITY ASSURANCE PROCEDURE

TS4
4.13.2011

QA SPECIALIST Cheryl STRAUB AGENCY Tucson Police Department
DATE 04/07/2011 TIME 05:59
INTOXILYZER SERIAL # 80-002759 LOCATION TEAM 4

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 097 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Cheryl Straub

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

297
CS

184

INTOXILYZER 8000
Location:
Serial Number : 80-002759
04/07/2011 06:01:15

Standard Lot#: 719919
Last Changed By: C. CONLEY #34517

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP,SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	06:02:10
Diagnostic Test	Pass	06:02:37
Air Blank	0.000	06:03:05
0.100 Cal Check	0.099	06:03:24
Air Blank	0.000	06:03:53
Subject Test	0.000	06:04:50
Air Blank	ABT*	06:04:57
Air Blank	0.000	06:05:25

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

377
CS

TS4

INTOXILYZER 8000
Location:
Serial Number : 80-002759
04/07/2011 06:05:29

Standard Lot#: 719919
Last Changed By: C. CONLEY #34517

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	06:06:01
Air Blank	IPS**	06:06:09

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

497
08

INTOXILYZER 8000
Location:
Serial Number : 80-002759
04/07/2011 06:06:18

BM

Standard Lot#: 719919
Last Changed By: C. CONLEY #34517

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	06:07:09
Diagnostic Test	Pass	06:07:36
Air Blank	0.000	06:08:04
0.100 Cal Check	0.098	06:08:22
Air Blank	0.000	06:08:52
Subject Test	0.071*	06:12:00
Air Blank	0.000	06:12:32
Air Blank	0.000	06:13:00
Subject Test	XXX**	06:13:47
Air Blank	0.000	06:14:21

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

577
c8

TS4

INTOXILYZER 8000
Location:
Serial Number: 80-002759
04/07/2011 06:14:43

Standard Lot#: 719919
Last Changed By: C. CONLEY #34517

QAS: C. STRAUB #50772
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	06:15:10
Cal Check	0.097	06:15:29
Air Blank	0.000	06:15:59



Operator's Signature

697
CS

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TS4
4.13.2011

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Cheryl Straub AGENCY Tucson Police Department
DATE 04/07/2011 TIME 05:59
Intoxilyzer Serial # 80-002759 Location TEAM 4

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution N/A AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Cheryl Straub

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

7g7
cs


154

INTOXILYZER 8000
Location:
Serial Number: 80-002759
04/07/2011 05:59:09

Standard Lot#: 719919
Last Changed By: C. CONLEY #34517

QAS: C. STRAUB #50772
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	05:59:36
Cal Check	0.099	05:59:55
Air Blank	0.000	06:00:25

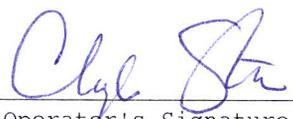

Operator's Signature