



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002759 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TSG
8.2.2011

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 07/28/11 Cal 1 0.095 AC Cal 2 0.094 AC

Between the dates of 06/30/11 and 07/28/11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Cheryl K. Straub

Cheryl K. Straub


08/01/11

ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 1 day of August, 2011.

Sue A. Soza

Notary Public

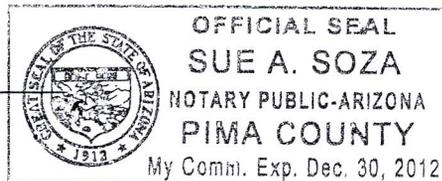

*Includes documents specified in A.A.C. R13-10-104.A.5

177
CS

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

TS4
8-2-2011

QA SPECIALIST Cheryl STRAUB AGENCY Tucson Police Department

DATE 07/28/11 TIME 05:23

INTOXILYZER SERIAL # 80-002759 LOCATION TEAM 4

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 094 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Cheryl Straub

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

277
c2

INTOXILYZER 8000

Location:

Serial Number : 80-002759

07/28/2011

05:25:55

TS4

Standard Lot#: 719919

Last Changed By: C. CONLEY

#34517

QAS: C. STRAUB

#50772

TUCSON PD

Operator: C. STRAUB

#50772

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	05:26:49
Diagnostic Test	Pass	05:27:16
Air Blank	0.000	05:27:44
0.100 Cal Check	0.096	05:28:02
Air Blank	0.000	05:28:32
Subject Test	0.000	05:28:52
Air Blank	ABT*	05:28:59
Air Blank	0.000	05:29:27

*Sequence Aborted

Not a Successfully Completed Test Sequence

377
CS

INTOXILYZER 8000

Location:

Serial Number : 80-002759

07/28/2011

05:29:31

TS4

Standard Lot#: 719919

Last Changed By: C. CONLEY

#34517

QAS: C. STRAUB

#50772

TUCSON PD

Operator: C. STRAUB

#50772

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	05:30:03
Air Blank	IPS**	05:30:11

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

477
05

INTOXILYZER 8000

Location:

Serial Number : 80-002759

07/28/2011

05:30:17

BY

Standard Lot#: 719919

Last Changed By: C. CONLEY

#34517

QAS: C. STRAUB

#50772

TUCSON PD

Operator: C. STRAUB

#50772

TUCSON PD

Subject: SQAP,SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	05:31:10
Diagnostic Test	Pass	05:31:37
Air Blank	0.000	05:32:04
0.100 Cal Check	0.097	05:32:23
Air Blank	0.000	05:32:52
Subject Test	0.400*	05:36:01
Air Blank	0.000	05:36:35
Air Blank	0.000	05:37:02
Subject Test	XXX**	05:37:39
Air Blank	0.000	05:38:13

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

677
09

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TS4
8.2.2011

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Cheryl STRAUB AGENCY Tucson Police Department
DATE 07/28/11 TIME 05:23
Intoxilyzer Serial # 80-002759 Location TEAM 4

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution N/A AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 095 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Cheryl Straub

1708



TUCSON POLICE DEPARTMENT CRIME LABORATORY

INTOXILYZER SERVICE RECORD FOR CALIBRATION GAS TANK CHANGE

709
8-2-2011

Date: 07/28/2011

QAS: Cheryl STRAUB

Intoxilyzer Serial Number: 80-002759

New Lot Number: 04911100A1

New Expiration Date: 03/01/2013

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."
Select "D" for dry gas and enter the following information:
 - Enter Target Value: 0.100
 - Enter Control Lot #: From the tank label
 - Enter Exp Date: From the tank label (mm/dd/yyyy)
 - Edit info Y/N? Y or N
 - Update Std Info? Y ENTER
 - "Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "Control Testing."
Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 098 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Cheryl Straub

272
CS

INTOXILYZER 8000

Location:

Serial Number: 80-002759

07/28/2011

05:50:10

TSK


Standard Lot#: 04911100A1

Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	05:50:37
Cal Check	0.098	05:50:56
Air Blank	0.000	05:51:25

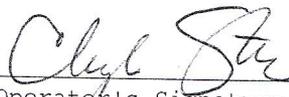
Operator's Signature