



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

SC
1/19/12

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002759 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 1.12.2012 Cal 1 0.098 AC Cal 2 0.096 AC

Between the dates of 12.16.2011 and 1.12.2012 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on .
The instrument was returned to service on .

Terry S Gallegos

Terry S. Gallegos


1.18.2012

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 18, day of January, 2012.

A.M. Hernandez

Notary Public

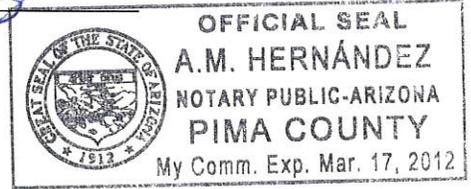

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7
TS4
SCV
1/19/12

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 1.12.2012 TIME 0947
INTOXILYZER SERIAL # 80-002759 LOCATION ODE

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0.000 AC.
(✓) 2. Error recognition logic system functioning. ✓✓✓
Not a Successfully Completed Test Sequence printed
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed. ✓
(✓) 4. Standard Calibration Check standard 0.100 AC
Result: 0.096 AC

Instrument is operating properly and accurately. YES ✓ NO

COMMENTS

SIGNATURE Terry Gallegos

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location:

Serial Number : 80-002759

01/12/2012

09:49:29

2077
TS4

Standard Lot#: 04911100A1

Last Changed By: C. STRAUB

#50772

SC

QAS: T. GALLEGOS

#28847

TUCSON PD

Operator: T. GALLEGOS

#28847

TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:50:23
Diagnostic Test	Pass	09:50:50
Air Blank	0.000	09:51:17
0.100 Cal Check	0.098	09:51:36
Air Blank	0.000	09:52:06
Subject Test	0.000	09:52:22
Air Blank	0.000	09:52:50
Five Minute Wait	ABT*	
Air Blank	0.000	09:53:23

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location:

Serial Number : 80-002759

01/12/2012

09:53:40

3077
754

Standard Lot#: 04911100A1

Last Changed By: C. STRAUB

#50772

QAS: T. GALLEGOS

#28847

TUCSON PD

Operator: T. GALLEGOS

#28847

TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:54:33
Diagnostic Test	Pass	09:55:00
Air Blank	RFI*	09:55:15
Air Blank	IPS**	09:55:23

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location:

Serial Number : 80-002759

01/12/2012

09:55:43

4 of 7
B4

Standard Lot#: 04911100A1

Last Changed By: C. STRAUB

#50772

QAS: T. GALLEGOS

#28847

TUCSON PD

5C

Operator: T. GALLEGOS

#28847

TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:56:37
Diagnostic Test	Pass	09:57:05
Air Blank	0.000	09:57:32
0.100 Cal Check	0.098	09:57:51
Air Blank	0.000	09:58:20
Subject Test	0.068*	10:01:29
Air Blank	0.000	10:02:00
Air Blank	0.000	10:02:27
Subject Test	XXX**	10:02:51
Air Blank	0.000	10:03:24

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location:

Serial Number: 80-002759

01/12/2012

10:03:48

5077

754

SC


Standard Lot#: 04911100A1

Last Changed By: C. STRAUB #50772

QAS: T. GALLEGOS #28847

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:04:15
Cal Check	0.096	10:04:34
Air Blank	0.000	10:05:03

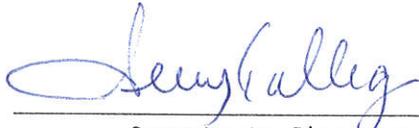

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

6077
TS4
SC ✓
1/19/12

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 1-12-2012 TIME 0947
Intoxilyzer Serial # 80-002759 Location ONE

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location:

Serial Number: 80-002759

01/12/2012

09:47:17

*7off
B4*

Standard Lot#: 04911100A1


Last Changed By: C. STRAUB #50772

SC

QAS: T. GALLEGOS #28847

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:47:44
Cal Check	0.098	09:48:03
Air Blank	0.000	09:48:33

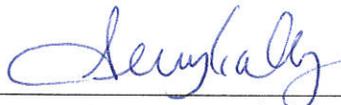

Operator's Signature