



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002759 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3-8-2012 Cal 1 0.098 AC Cal 2 0.098 AC

Between the dates of 2-9-2012 and 3-8-2012 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -
The instrument was returned to service on -

[Handwritten signature of Terry S. Gallegos]

Terry S. Gallegos

3-26-2012

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 26th day of March, 2012.

[Handwritten signature of Elma Escarcega]

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

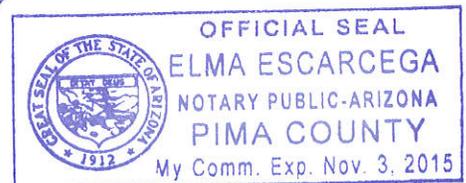


EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7
TS4
SC
3/27/12

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 3.8.2012 TIME 0857

INTOXILYZER SERIAL # 80-002759 LOCATION ODE

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 098 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Terry Gallegos

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location:

Serial Number : 80-002759

03/08/2012

08:59:03

2 of 7
TS4

Standard Lot#: 04911100A1

Last Changed By: C. STRAUB

#50772

5C

QAS: T. GALLEGOS

#28847

TUCSON PD

Operator: T. GALLEGOS

#28847

TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:59:57
Diagnostic Test	Pass	09:00:25
Air Blank	0.000	09:00:52
0.100 Cal Check	0.097	09:01:11
Air Blank	0.000	09:01:40
Subject Test	0.000	09:02:04
Air Blank	ABT*	09:02:15
Air Blank	0.000	09:02:43

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location:

Serial Number : 80-002759

03/08/2012

09:03:09

30f7
TJ4

Standard Lot#: 04911100A1

Last Changed By: C. STRAUB

#50772

SC

QAS: T. GALLEGOS

#28847

TUCSON PD

Operator: T. GALLEGOS

#28847

TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:04:02
Diagnostic Test	Pass	09:04:29
Air Blank	RFI*	09:04:44
Air Blank	IPS**	09:04:53

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location:

Serial Number : 80-002759

03/08/2012

09:05:11

4087

TJ4

3C

Standard Lot#: 04911100A1

Last Changed By: C. STRAUB

#50772

QAS: T. GALLEGOS

#28847

TUCSON PD

Operator: T. GALLEGOS

#28847

TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:06:02
Diagnostic Test	Pass	09:06:29
Air Blank	0.000	09:06:56
0.100 Cal Check	0.098	09:07:15
Air Blank	0.000	09:07:45
Subject Test	0.225*	09:10:53
Air Blank	0.000	09:11:25
Air Blank	0.000	09:11:52
Subject Test	XXX**	09:12:11
Air Blank	0.000	09:12:45

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6087
TJ4
SC
3/27/12

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 3-8-2012 TIME 0857
Intoxilyzer Serial # 80-002759 Location ODE

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

