



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002759 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

SC  
7/9/12

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6.28.2012 Cal 1 0.097 AC Cal 2 0.097 AC

Between the dates of 5.31.2012 and 6.28.2012 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
The instrument was returned to service on -.

Terry S. Gallegos

Terry S. Gallegos


7.3.2012

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 3<sup>rd</sup> day of July, 2012.

Adina Escarcega

Notary Public

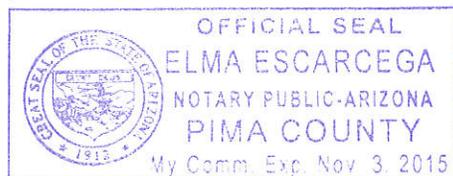

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7  
TS4  
SC  
7/9/12

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 6-28-2012 TIME 0920

INTOXILYZER SERIAL # 80-002759 LOCATION ODE

(  ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (  ) 1. Clock time check.
- (  ) 2. Date check.

OPERATIONAL TESTS

- (  ) 1. Alcohol-free subject test result 0. 000 AC.
- (  ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (  ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (  ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 097 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE 

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location:

Serial Number : 80-002759

06/28/2012

09:26:06

20f7  
TS4

---

Standard Lot#: 23511100A2

Last Changed By: T. GALLEGOS

#28847

SC

---

QAS: T. GALLEGOS

#28847

TUCSON PD

---

Operator: T. GALLEGOS

#28847

TUCSON PD

---

Subject: TEST, SQAP

DOB: 11/11/1987

Sex: F

Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	09:27:03
Diagnostic Test	Pass	09:27:30
Air Blank	0.000	09:27:58
0.100 Cal Check	0.098	09:28:17
Air Blank	0.000	09:28:46
Subject Test	0.000	09:29:17
Air Blank	0.000	09:29:44
Five Minute Wait	ABT*	
Air Blank	0.000	09:30:36

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000

Location:

Serial Number : 80-002759

06/28/2012

09:30:43

3 of 7  
TS4

Standard Lot#: 23511100A2

Last Changed By: T. GALLEGOS

#28847

SC

QAS: T. GALLEGOS

#28847

TUCSON PD

Operator: T. GALLEGOS

#28847

TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1987

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	09:31:36
Air Blank	IPS**	09:31:45

\*RFI Detect

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location:  
Serial Number : 80-002759  
06/28/2012 09:32:01

40F7  
TSh

Standard Lot#: 23511100A2  
Last Changed By: T. GALLEGOS #28847

5C

QAS: T. GALLEGOS #28847  
TUCSON PD

Operator: T. GALLEGOS #28847  
TUCSON PD

Subject: TEST,SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:32:50
Diagnostic Test	Pass	09:33:17
Air Blank	0.000	09:33:45
0.100 Cal Check	0.098	09:34:03
Air Blank	0.000	09:34:33
Subject Test	0.247*	09:37:41
Air Blank	0.000	09:38:14
Air Blank	0.000	09:38:41
Subject Test	XXX**	09:39:01
Air Blank	0.000	09:39:35

\*Deficient Sample  
\*\*Invalid Sample  
Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000

Location:

Serial Number: 80-002759

06/28/2012

09:42:48

50F7  
TS4

-----  
Standard Lot#: 23511100A2

Last Changed By: T. GALLEGOS #28847


SC

-----  
QAS: T. GALLEGOS #28847

TUCSON PD

-----  
-----

Test	g/210L	Time
Air Blank	0.000	09:43:16
Cal Check	0.097	09:43:34
Air Blank	0.000	09:44:04

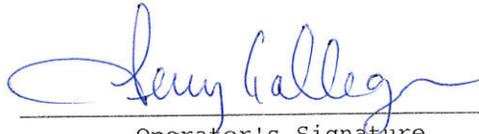

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000


6077  
TS6

SC  
7/9/12

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 6-28-2012 TIME 0920  
Intoxilyzer Serial # 80-002759 Location ODE

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 097 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

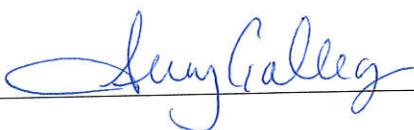SIGNATURE 

DPS Form Exh G-2 (Rev. 05-01)

