



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

SC  
7/31/12

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002759 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 7.26.2012 Cal 1 0.098 AC Cal 2 0.098 AC

Between the dates of 6.28.2012 and 7.26.2012 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.  
The instrument was returned to service on —.

Terry Gallegos

Terry S. Gallegos


7.31.2012

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 31, day of July, 2012.

Ann L. Bercroft

Notary Public

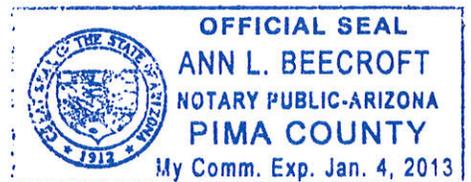

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

10 of 7

TS4

SC  
7/31/12

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 7.26.2012 TIME 0906

INTOXILYZER SERIAL # 80-002759 LOCATION ODE

(  ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (  ) 1. Clock time check.
- (  ) 2. Date check.

OPERATIONAL TESTS

- (  ) 1. Alcohol-free subject test result 0.000 AC.
- (  ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (  ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (  ) 4. Standard Calibration Check standard 0.100 AC  
Result: 0.098 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE Terry Gallegos

INTOXILYZER 8000

Location:

Serial Number : 80-002759

07/26/2012

09:10:01

Standard Lot#: 23511100A2

Last Changed By: T. GALLEGOS

#28847

QAS: T. GALLEGOS

#28847

TUCSON PD

Operator: T. GALLEGOS

#28847

TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1987

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:10:58
Diagnostic Test	Pass	09:11:25
Air Blank	0.000	09:11:52
0.100 Cal Check	0.099	09:12:11
Air Blank	0.000	09:12:40
Subject Test	0.000	09:13:09
Air Blank	0.000	09:13:37
Five Minute Wait	ABT*	
Air Blank	0.000	09:14:12

\*Sequence Aborted

2077  
TS4

SC

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000

Location:

Serial Number : 80-002759

07/26/2012

09:24:31

*TS4*  
*307*      *SC*

Standard Lot#: 23511100A2

Last Changed By: T. GALLEGOS

#28847

QAS: T. GALLEGOS

#28847

TUCSON PD

Operator: T. GALLEGOS

#28847

TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1987

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:25:22
Diagnostic Test	Fail*	09:25:28
Air Blank	IPS**	09:25:42

\* Analytical Stability Test

\*\*Improper Sample

*→ RFI during diagnostics  
Provided by officer Gamble  
(PR 44527)*

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000

Location:

Serial Number : 80-002759

07/26/2012

09:14:24

4077  
Bh

Standard Lot#: 23511100A2

Last Changed By: T. GALLEGOS

#28847

SC

QAS: T. GALLEGOS

#28847

TUCSON PD

Operator: T. GALLEGOS

#28847

TUCSON PD

Subject: TEST,SQAP

DOB: 11/11/1987

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:15:33
Diagnostic Test	Pass	09:16:00
Air Blank	0.000	09:16:28
0.100 Cal Check	0.098	09:16:46
Air Blank	0.000	09:17:16
Subject Test	0.232*	09:20:24
Air Blank	0.000	09:20:56
Air Blank	0.000	09:21:24
Subject Test	XXX**	09:21:43
Air Blank	0.000	09:22:16

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**



EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6077  
TS4  
SC  
7/31/12

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 7.26.2012 TIME 0906  
Intoxilyzer Serial # 80-002759 Location ODE

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 098 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location:

Serial Number: 80-002759

07/26/2012

09:06:59

70f7  
TS4

-----  
Standard Lot#: 23511100A2


Last Changed By: T. GALLEGOS #28847

-----  
QAS: T. GALLEGOS #28847

TUCSON PD

SC

Test	g/210L	Time
Air Blank	0.000	09:07:27
Cal Check	0.098	09:07:45
Air Blank	0.000	09:08:15

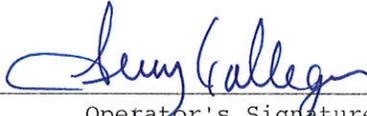

Operator's Signature