



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

SC  
3/1/13

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002759 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 2.7.2013 Cal 1 0.097 AC Cal 2 0.096 AC

Between the dates of 1.9.2013 and 2.7.2013 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
The instrument was returned to service on -.

Terry S. Gallegos

Terry S. Gallegos

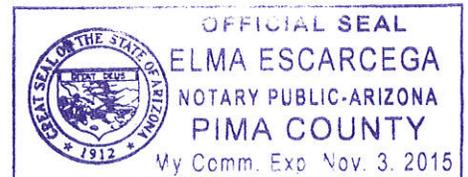
2.26.2013

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 26<sup>th</sup>, day of February, 2013.

Elma Escarcega

Notary Public



\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1 of 7  
TS4

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD QUALITY ASSURANCE PROCEDURE

SC  
3/1/13

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 2-7-2013 TIME 0856  
INTOXILYZER SERIAL # 80-002759 LOCATION ODE

(  ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (  ) 1. Clock time check.
- (  ) 2. Date check.

OPERATIONAL TESTS

- (  ) 1. Alcohol-free subject test result 0. 000 AC.
- (  ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (  ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (  ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 096 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE

Terry Gallegos

INTOXILYZER 8000

Location:

Serial Number : 80-002759

02/07/2013

08:57:48

2 of 7  
TSG

Standard Lot#: 23511100A2

Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847

TUCSON PD

Operator: T. GALLEGOS #28847

TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1987

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:58:42
Diagnostic Test	Pass	08:59:09
Air Blank	0.000	08:59:37
0.100 Cal Check	0.097	08:59:55
Air Blank	0.000	09:00:25
Subject Test	0.000	09:00:48
Air Blank	0.000	09:01:16
Five Minute Wait	ABT*	
Air Blank	0.000	09:01:54

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

30f7  
TSh

INTOXILYZER 8000  
Location:  
Serial Number : 80-002759  
02/07/2013 09:02:00

5c

Standard Lot#: 23511100A2  
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847  
TUCSON PD

Operator: T. GALLEGOS #28847  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:02:53
Diagnostic Test	Pass	09:03:21
Air Blank	0.000	09:03:48
0.100 Cal Check	0.097	09:04:07
Air Blank	RFI*	09:04:14
Air Blank	IPS**	09:04:23

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

4 of 7  
DJ

INTOXILYZER 8000  
Location:  
Serial Number : 80-002759  
02/07/2013

09:04:32

5<

Standard Lot#: 23511100A2  
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847  
TUCSON PD

Operator: T. GALLEGOS #28847  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:05:23
Diagnostic Test	Pass	09:05:50
Air Blank	0.000	09:06:18
0.100 Cal Check	0.097	09:06:36
Air Blank	0.000	09:07:06
Subject Test	0.186*	09:10:14
Air Blank	0.000	09:10:47
Air Blank	0.000	09:11:14
Subject Test	XXX**	09:11:36
Air Blank	0.000	09:12:10

\*Deficient Sample  
\*\*Invalid Sample  
Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000

Location:

Serial Number: 80-002759

02/07/2013

09:12:30

5087  
By

5c

-----  
Standard Lot#: 23511100A2

Last Changed By: T. GALLEGOS #28847

-----  
QAS: T. GALLEGOS #28847

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:12:58
Cal Check	0.096	09:13:16
Air Blank	0.000	09:13:46

  
-----  
Operator's Signature

EXHIBIT G-2

6077  
TS4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

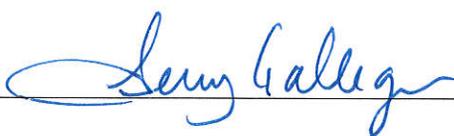
SC  
3/1/13

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 2-7-2013 TIME 0856  
Intoxilyzer Serial # 80-002759 Location ODE

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 097 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE 

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location:

Serial Number: 80-002759

02/07/2013

08:56:23

70f7  
TS4

5c

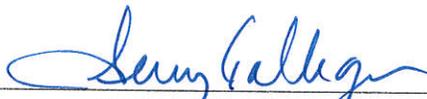
-----  
Standard Lot#: 23511100A2

Last Changed By: T. GALLEGOS #28847

-----  
QAS: T. GALLEGOS #28847

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	08:56:50
Cal Check	0.097	08:57:09
Air Blank	0.000	08:57:38

  
-----  
Operator's Signature