



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

SC  
3/22/13

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002759 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3-7-2013 Cal 1 0.098 AC Cal 2 0.097 AC

Between the dates of 2-7-2013 and 3-7-2013 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.  
The instrument was returned to service on —.


Terry S. Gallegos

3-13-2013

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 13<sup>th</sup>, day of March, 2013.

Notary Public

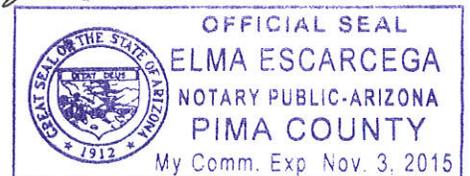

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1087  
T24

SC  
3/22/13

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 3-7-2013 TIME 0921

INTOXILYZER SERIAL # 80-002759 LOCATION ODE

(  ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (  ) 1. Clock time check.
- (  ) 2. Date check.

OPERATIONAL TESTS

- (  ) 1. Alcohol-free subject test result 0. 000 AC.
- (  ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (  ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (  ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 097 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Terry Gallegos

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location:

Serial Number : 80-002759

03/07/2013

09:23:14

2 of 7  
TB4

Standard Lot#: 23511100A2

Last Changed By: T. GALLEGOS #28847

SC

QAS: T. GALLEGOS #28847

TUCSON PD

Operator: T. GALLEGOS #28847

TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1987

Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:24:07
Diagnostic Test	Pass	09:24:34
Air Blank	0.000	09:25:02
0.100 Cal Check	0.098	09:25:20
Air Blank	0.000	09:25:50
Subject Test	0.000	09:26:44
Air Blank	0.000	09:27:11
Five Minute Wait	ABT*	
Air Blank	0.000	09:27:50

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000

Location:

Serial Number : 80-002759

03/07/2013

09:27:56

30f7  
TS4

Standard Lot#: 23511100A2

Last Changed By: T. GALLEGOS #28847

SC

QAS: T. GALLEGOS #28847

TUCSON PD

Operator: T. GALLEGOS #28847

TUCSON PD

Subject: TEST,SQAP

DOB: 11/11/1987

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:28:48
Diagnostic Test	Pass	09:29:15
Air Blank	0.000	09:29:43
0.100 Cal Check	0.098	09:30:02
Air Blank	RFI*	09:30:10
Air Blank	IPS**	09:30:20

\*RFI Detect

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000

Location:

Serial Number : 80-002759

03/07/2013

09:30:39

4 of 7  
TS4

Standard Lot#: 23511100A2

Last Changed By: T. GALLEGOS

#28847

5<

QAS: T. GALLEGOS

#28847

TUCSON PD

Operator: T. GALLEGOS

#28847

TUCSON PD

Subject: TEST,SOAP

DOB: 11/11/1987

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:31:31
Diagnostic Test	Pass	09:31:59
Air Blank	0.000	09:32:26
0.100 Cal Check	0.099	09:32:45
Air Blank	0.000	09:33:14
Subject Test	0.117*	09:36:23
Air Blank	0.000	09:36:54
Air Blank	0.000	09:37:22
Subject Test	XXX**	09:37:38
Air Blank	0.000	09:38:12

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**



EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

6 of 7  
TJG

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

SC  
3/22/13

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 3.7.2013 TIME 0921  
Intoxilyzer Serial # 80-002759 Location ODE

- (  ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (  ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (  ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (  ) 4. Air blank completed.
- (  ) 5. Calibration check completed. Test results 0. 098 AC.
- (  ) 6. Air blank completed.
- (  ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)

