



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

SC  
2/14/14

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002759 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 2.6.2014 Cal 1 0.096 AC Cal 2 0.096 AC

Between the dates of 1.9.2014 and 2.6.2014 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.  
The instrument was returned to service on —.

Terry S. Gallegos

Terry S. Gallegos

2.13.2014

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 13<sup>th</sup>, day of February, 2014.

Elma Escarcega

Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

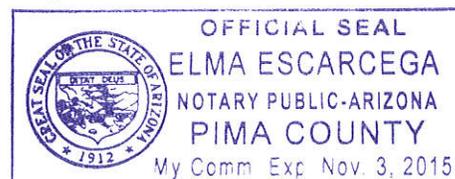


EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1057  
Dy SC  
2/14/14

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 2-6-2014 TIME 0920  
INTOXILYZER SERIAL # 80-002759 LOCATION ODE

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 096 AC

Instrument is operating properly and accurately. YES ✓ NO \_\_\_\_\_

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Terry Gallegos

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002759  
Core Version: 8105.48  
02/06/2014

09:22:04

*dof 7*  
*TS4*

Standard Lot#: 30412100A1  
Last Changed By: T. GALLEGOS #28847

*SL*

QAS: T. GALLEGOS #28847  
TUCSON PD

Operator: T. GALLEGOS #28847  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:22:53
Diagnostic Test	Pass	09:23:28
Air Blank	0.000	09:23:56
0.100 Cal Check	0.097	09:24:16
Air Blank	0.000	09:24:46
Subject Test	0.000	09:25:14
Air Blank	0.000	09:25:42
Wait	ABT*	
Air Blank	0.000	09:26:20

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002759  
Core Version: 8105.48  
02/06/2014

09:26:52

3057  
TS4

5C

Standard Lot#: 30412100A1  
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847  
TUCSON PD

Operator: T. GALLEGOS #28847  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:27:43
Diagnostic Test	Pass	09:28:18
Air Blank	0.000	09:28:46
0.100 Cal Check	RFI*	09:29:05
Air Blank	IPS**	09:29:13

\*RFI Detect

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002759  
Core Version: 8105.48  
02/06/2014

09:29:20

4 of 7  
By

Standard Lot#: 30412100A1  
Last Changed By: T. GALLEGOS #28847

50

QAS: T. GALLEGOS #28847  
TUCSON PD

Operator: T. GALLEGOS #28847  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:30:10
Diagnostic Test	Pass	09:30:45
Air Blank	0.000	09:31:12
0.100 Cal Check	0.097	09:31:33
Air Blank	0.000	09:32:02
Subject Test	0.184*	09:35:10
Air Blank	0.000	09:35:42
Air Blank	0.000	09:36:09
Subject Test	XXX**	09:36:37
Air Blank	0.000	09:37:11

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002759  
Core Version: 8105.48  
02/06/2014 09:37:29

507  
TS4  
5C

-----  
Standard Lot#: 30412100A1  
Last Changed By: T. GALLEGOS #28847

-----  
QAS: T. GALLEGOS #28847  
TUCSON PD

-----  
Operator: T. GALLEGOS #28847  
TUCSON PD

-----  
Test g/210L Time  
-----  
Air Blank 0.000 09:37:59  
0.100 Cal Check 0.096 09:38:20  
Air Blank 0.000 09:38:50

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7  
TSG  
SC  
2/14/14

QA Specialist Terry Gallegos 20887 AGENCY Tucson Police Department  
DATE 2-6-2014 TIME 0920  
Intoxilyzer Serial # 80-002759 Location ODE

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 096 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE 

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002759  
Core Version: 8105.48  
02/06/2014 09:20:26

7077  
154

-----  
Standard Lot#: 30412100A1  
Last Changed By: T. GALLEGOS #28847  
-----

5<

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB  
-----

Test	g/210L	Time
Air Blank	0.000	09:20:53
0.100 Cal Check	0.096	09:21:14
Air Blank	0.000	09:21:43