



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002759 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

SL
3/20/14

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3.6.2014 Cal 1 0.096 AC Cal 2 0.096 AC

Between the dates of 2.6.2014 and 3.6.2014 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -
The instrument was returned to service on -

[Handwritten signature of Terry S. Gallegos]

Terry S. Gallegos


3.19.2014

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 19th, day of March, 2014.

[Handwritten signature of Notary Public]

Notary Public

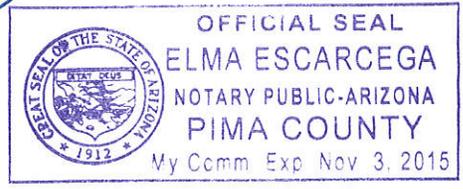

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

10f7
TS4

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

SC
3/20/14

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 3.6.2014 TIME 0919

INTOXILYZER SERIAL # 80-002759 LOCATION ODE

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 096 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE 

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002759
Core Version: 8105.48
03/06/2014

09:21:39

2 of 7
TS4
SL

Standard Lot#: 30412100A1
Last Changed By: T. GALLEGOS #28847

OAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SOAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:22:32
Diagnostic Test	Pass	09:23:07
Air Blank	0.000	09:23:34
0.100 Cal Check	0.096	09:23:55
Air Blank	0.000	09:24:24
Subject Test	0.000	09:24:56
Air Blank	0.000	09:25:23
Wait	ABT*	
Air Blank	0.000	09:25:59

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002759
Core Version: 8105.48
03/06/2014

09:26:48

30 of 7
TSh
SL

Standard Lot#: 30412100A1
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:27:43
Diagnostic Test	Pass	09:28:18
Air Blank	0.000	09:28:46
0.100 Cal Check	0.097	09:29:07
Air Blank	RFI*	09:29:15
Air Blank	IPS**	09:29:25

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002759
Core Version: 8105.48
03/06/2014

09:29:30

4 of 7
154

SC

Standard Lot#: 30412100A1
Last Changed By: T. GALLEGOS #28847

OAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:30:19
Diagnostic Test	Pass	09:30:54
Air Blank	0.000	09:31:21
0.100 Cal Check	0.097	09:31:42
Air Blank	0.000	09:32:11
Subject Test	0.124*	09:35:19
Air Blank	0.000	09:35:51
Air Blank	0.000	09:36:18
Subject Test	XXX**	09:36:54
Air Blank	0.000	09:37:27

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002759
Core Version: 8105.48
03/06/2014 09:39:49

50f 7
T34

5C

Standard Lot#: 30412100A1
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:40:19
0.100 Cal Check	0.096	09:40:40
Air Blank	0.000	09:41:09

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6077
B4

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

5C
3/20/14

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 3.6.2014 TIME 0919
Intoxilyzer Serial # 80-002759 Location ODE

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 096 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE 

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002759
Core Version: 8105.48
03/06/2014 09:19:18

70f7
TS4

Standard Lot#: 30412100A1
Last Changed By: T. GALLEGOS #28847

5L

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	09:19:46
0.100 Cal Check	0.096	09:20:06
Air Blank	0.000	09:20:36