



CERTIFICATION by CUSTODIAN of RECORDS  
 Tucson Police Department Crime Laboratory  
 1306 W. Miracle Mile, Tucson, AZ 85705

SC  
4/22/14

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002759 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 4.3.2014 Cal 1 0.096 AC Cal 2 0.097 AC

Between the dates of 3.6.2014 and 4.3.2014 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.  
 The instrument was returned to service on —.

Terry S. Gallegos

Terry S. Gallegos


4.22.2014

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 22nd, day of April, 2014.

Edma Escarcega

Notary Public

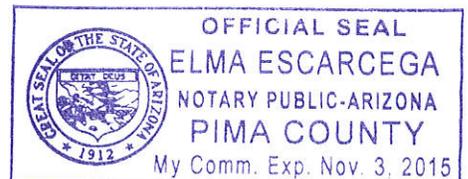

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

10f9  
TSK  
SC  
4/22/14

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 4.3.2014 TIME 0903

INTOXILYZER SERIAL # 80-002759 LOCATION ODE

(  ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (  ) 1. Clock time check.
- (  ) 2. Date check.

OPERATIONAL TESTS

- (  ) 1. Alcohol-free subject test result 0. 000 AC.
- (  ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (  ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (  ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 097 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE 

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002759  
Core Version: 8105.48  
04/03/2014

09:04:57

2 of 9  
BH  
SC

---

Standard Lot#: 30412100A1  
Last Changed By: T. GALLEGOS #28847

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 09:06:00 |
| Diagnostic Test | Pass   | 09:06:35 |
| Air Blank       | 0.000  | 09:07:03 |
| 0.100 Cal Check | 0.097  | 09:07:23 |
| Air Blank       | 0.000  | 09:07:53 |
| Subject Test    | 0.000  | 09:08:25 |
| Air Blank       | 0.000  | 09:08:53 |
| Wait            | ABT*   |          |
| Air Blank       | 0.000  | 09:09:27 |

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002759  
Core Version: 8105.48  
04/03/2014

09:09:37

3 of 9  
D4  
SC

---

Standard Lot#: 30412100A1  
Last Changed By: T. GALLEGOS #28847

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST,SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 09:10:34 |
| Diagnostic Test | Pass   | 09:11:10 |
| Air Blank       | 0.000  | 09:11:37 |
| 0.100 Cal Check | RFI*   | 09:11:53 |
| Air Blank       | IPS**  | 09:12:02 |

\*RFI Detect

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002759  
Core Version: 8105.48  
04/03/2014

09:12:54

4 of 9  
D4

5C

Standard Lot#: 30412100A1  
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847  
TUCSON PD

Operator: T. GALLEGOS #28847  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 09:13:43 |
| Diagnostic Test | Pass   | 09:14:18 |
| Air Blank       | 0.000  | 09:14:46 |
| 0.100 Cal Check | 0.097  | 09:15:07 |
| Air Blank       | 0.000  | 09:15:36 |
| Subject Test    | 0.000* | 09:18:43 |
| Air Blank       | 0.000  | 09:19:11 |
| Air Blank       | 0.000  | 09:19:38 |
| Subject Test    | 0.000* | 09:22:46 |
| Air Blank       | 0.000  | 09:23:14 |

\*Deficient Sample

*Did not have breath spray with alcohol for these tests.*  
TS4

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002759  
Core Version: 8105.48  
04/03/2014

09:26:54

50f9  
TB4  
SC

Standard Lot#: 30412100A1  
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847  
TUCSON PD

Operator: T. GALLEGOS #28847  
TUCSON PD

Subject: TEST,SOAP  
DOB: 11/11/1987  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 09:27:47 |
| Diagnostic Test | Pass   | 09:28:22 |
| Air Blank       | 0.000  | 09:28:49 |
| 0.100 Cal Check | 0.097  | 09:29:10 |
| Air Blank       | 0.000  | 09:29:39 |
| Subject Test    | RGE*   | 09:29:59 |
| Air Blank       | 0.000  | 09:30:34 |

\*Range Exceeded

*Used hairspray on the mouthpiece for Invalid sample.  
TB4*

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002759  
Core Version: 8105.48  
04/03/2014

09:30:43

60f9  
TSh

5C

Standard Lot#: 30412100A1  
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847  
TUCSON PD

Operator: T. GALLEGOS #28847  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 09:31:35 |
| Diagnostic Test | Pass   | 09:32:10 |
| Air Blank       | 0.000  | 09:32:37 |
| 0.100 Cal Check | 0.097  | 09:32:58 |
| Air Blank       | 0.000  | 09:33:28 |
| Subject Test    | XXX*   | 09:33:48 |
| Air Blank       | 0.000  | 09:34:20 |

\*Invalid Sample  
Begin new deprivation period

Used hairspray with alcohol on the mouthpiece.  
TSh

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002759  
Core Version: 8105.48  
04/03/2014 09:34:36

70F9  
TS4

-----  
Standard Lot#: 30412100A1  
Last Changed By: T. GALLEGOS #28847

SC

-----  
QAS: T. GALLEGOS #28847  
TUCSON PD

-----  
Operator: T. GALLEGOS #28847  
TUCSON PD

-----  
Test g/210L Time  
-----  
Air Blank 0.000 09:35:06  
0.100 Cal Check 0.097 09:35:27  
Air Blank 0.000 09:35:57

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

80f9  
D4

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

5C  
4/22/14

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 4-3-2014 TIME 0903  
Intoxilyzer Serial # 80-002759 Location ODE

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 096 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE 

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002759  
Core Version: 8105.48  
04/03/2014 09:03:29

9089  
D4  
5c

-----  
Standard Lot#: 30412100A1  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB  
-----

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 09:03:56 |
| 0.100 Cal Check | 0.096  | 09:04:17 |
| Air Blank       | 0.000  | 09:04:47 |