



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

SC
6/6/14

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002759 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 5.29.2014 Cal 1 0.096 AC Cal 2 0.095 AC

Between the dates of 5.1.2014 and 5.29.2014 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.
The instrument was returned to service on —.

Terry S. Gallegos

Terry S. Gallegos


6.4.2014

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 24th, day of June, 2014.

Elma Escarcega

Notary Public

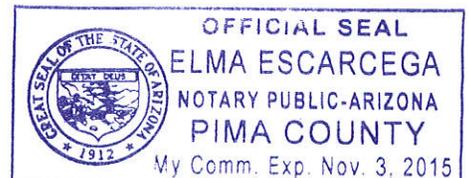

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1078
TS4
sc
6/2/14

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 5.29.2014 TIME 0914
INTOXILYZER SERIAL # 80-002759 LOCATION ODE

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 095 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE

Terry Gallegos

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002759
Core Version: 8105.48
05/29/2014 09:15:45

2078
TS4
SC

Standard Lot#: 30412100A1
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:16:36
Diagnostic Test	Pass	09:17:11
Air Blank	0.000	09:17:38
0.100 Cal Check	0.097	09:17:59
Air Blank	0.000	09:18:28
Subject Test	0.000	09:18:52
Air Blank	0.000	09:19:19
Wait	ABT*	
Air Blank	0.000	09:19:53

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002759
Core Version: 8105.48
05/29/2014

09:20:15

30f8
TS4
3C

Standard Lot#: 30412100A1
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TESTS, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:21:04
Diagnostic Test	Pass	09:21:39
Air Blank	0.000	09:22:06
0.100 Cal Check	RFI*	09:22:22
Air Blank	IPS**	09:22:31

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002759
Core Version: 8105.48
05/29/2014

09:22:51

40f8
TS4

3C

Standard Lot#: 30412100A1
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TESTS, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:23:44
Diagnostic Test	Pass	09:24:19
Air Blank	0.000	09:24:47
0.100 Cal Check	0.096	09:25:07
Air Blank	0.000	09:25:37
Subject Test	0.205*	09:28:44
Air Blank	0.000	09:29:16
Air Blank	0.000	09:29:44
Subject Test	0.259**	09:30:07
Air Blank	0.000	09:30:47
Wait		
Air Blank	0.000	09:35:34
Subject Test	0.000**	09:36:02
Air Blank	0.000	09:36:30

*Deficient Sample

**No .020 agreement

*Ran out of breath trying to do Invalid Sample. Waited 5 minutes,
then did another test. TS4*

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002759
Core Version: 8105.48
05/29/2014

09:36:34

50f8
TS4

5c

Standard Lot#: 30412100A1
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST,SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:37:24
Diagnostic Test	Pass	09:37:59
Air Blank	0.000	09:38:26
0.100 Cal Check	0.097	09:38:47
Air Blank	0.000	09:39:17
Subject Test	XXX*	09:39:44
Air Blank	0.000	09:40:18

*Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002759
Core Version: 8105.48
05/29/2014 09:40:41

60f8
TSG

3C

Standard Lot#: 30412100A1
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:41:10
0.100 Cal Check	0.095	09:41:31
Air Blank	0.000	09:42:00

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

70f8
TS4

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

SC
6/6/14

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 5.29.2014 TIME 0914
Intoxilyzer Serial # 80-002759 Location ODE

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0.100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0.096 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002759
Core Version: 8105.48
05/29/2014 09:14:07

80f8
TS4

Standard Lot#: 30412100A1
Last Changed By: T. GALLEGOS #28847

5<

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	09:14:35
0.100 Cal Check	0.096	09:14:55
Air Blank	0.000	09:15:25