



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TS4
1.21.2015

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002759 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 1/8/2015 Cal 1 0.099 AC Cal 2 0.098 AC

Between the dates of 12/15/2014 and 1/8/2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on _____.
The instrument was returned to service on _____.

William L Loomis

William L. Loomis


1/8/2015

ADPS QAS Permit # 33720

SUBSCRIBED AND SWORN to before me this 8th, day of January, 2015.

Elma Escarcega

Notary Public

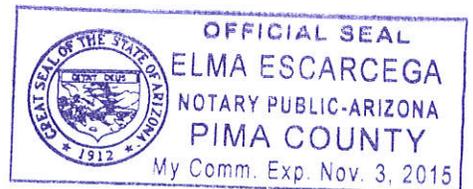

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7 WLL
By
1-21-2015

QA SPECIALIST Will Loomis 53545 AGENCY Tucson Police Department

DATE 1/8/2015 TIME 12:59

INTOXILYZER SERIAL # 80-002759 LOCATION CRIME LAB

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 098 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Will Loomis

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002759
Core Version: 8105.48
01/08/2015

13:01:07

2 of 7
WLL
TSY

Standard Lot#: 11514100A2
Last Changed By: S. CIONE #50616

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:02:09
Diagnostic Test	Pass	13:02:44
Air Blank	0.000	13:03:11
0.100 Cal Check	0.099	13:03:32
Air Blank	0.000	13:04:01
Subject Test	0.000	13:04:17
Air Blank	0.000	13:04:44
Wait	ABT*	
Air Blank	0.000	13:05:19

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002759
Core Version: 8105.48
01/08/2015

13:05:33

3 of 7
WLL
TS4

Standard Lot#: 11514100A2
Last Changed By: S. CIONE #50616

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:06:32
Diagnostic Test	Pass	13:07:07
Air Blank	0.000	13:07:35
0.100 Cal Check	RFI*	13:07:50
Air Blank	IPS**	13:07:59

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002759
Core Version: 8105.48
01/08/2015

13:08:07

4 of 7
WLL
Tsh

Standard Lot#: 11514100A2
Last Changed By: S. CIONE #50616

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:09:02
Diagnostic Test	Pass	13:09:37
Air Blank	0.000	13:10:04
0.100 Cal Check	0.099	13:10:25
Air Blank	0.000	13:10:54
Subject Test	0.131*	13:14:02
Air Blank	0.000	13:14:34
Air Blank	0.000	13:15:01
Subject Test	XXX**	13:15:25
Air Blank	0.000	13:15:59

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002759
Core Version: 8105.48
01/08/2015 13:16:27

5 of 7
WLL
T34

Standard Lot#: 11514100A2
Last Changed By: S. CIONE #50616

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	13:16:56
0.100 Cal Check	0.098	13:17:17
Air Blank	0.000	13:17:47

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7
will
TSK
1.21.2015

QA Specialist Will Loomis 53545 AGENCY Tucson Police Department
DATE 1/8/2015 TIME 12:59
Intoxilyzer Serial # 80-002759 Location CRIME LAB

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Will Loomis

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002759
Core Version: 8105.48
01/08/2015 12:59:36

7 of 7
WLL

Standard Lot#: 11514100A2
Last Changed By: S. CIONE #50616

TS4

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	13:00:07
0.100 Cal Check	0.099	13:00:27
Air Blank	0.000	13:00:57

Preventative Maintenance
Performed -- Timer Reset