



CERTIFICATION by CUSTODIAN of RECORDS  
 Tucson Police Department Crime Laboratory  
 1306 W. Miracle Mile, Tucson, AZ 85705

TS4  
 3-11-2015


I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002759 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3-5-15 Cal 1 0.099 AC Cal 2 0.098 AC

Between the dates of 2-6-15 and 3-5-15 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.  
 The instrument was returned to service on N/A.

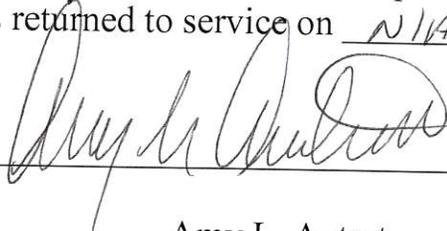  
 \_\_\_\_\_  
 Amy L. Autret


3-11-15

ADPS QAS Permit # 33718

SUBSCRIBED AND SWORN to before me this 11<sup>th</sup> day of March, 2015.

  
 \_\_\_\_\_

Notary Public

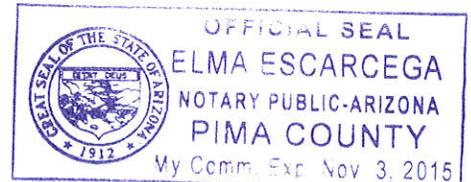

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1 of 7  
dB  
B4  
3.11.2015

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Amy Autret 53546 AGENCY Tucson Police Department

DATE 3/5/15 TIME 1306

INTOXILYZER SERIAL # 80-002759 LOCATION OPW

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS


- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 098 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE

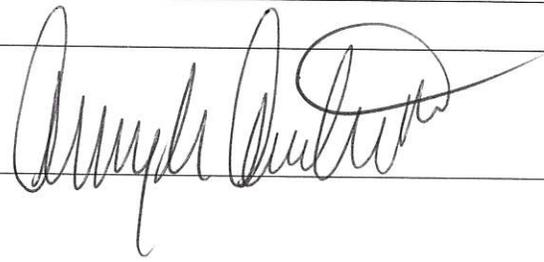

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2 of 7  
[Signature]

T34

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002759  
Core Version: 8105.48  
03/05/2015 13:10:09

Standard Lot#: 11514100A2  
Last Changed By: S. CIONE #50616

QAS: A. AUTRET #53546  
TUCSON PD

Operator: A. AUTRET #53546  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:11:30
Diagnostic Test	Pass	13:12:05
Air Blank	0.000	13:12:33
0.100 Cal Check	0.099	13:12:53
Air Blank	0.000	13:13:23
Subject Test	0.000	13:13:43
Air Blank	0.000	13:14:11
Wait	ABT*	
Air Blank	0.000	13:14:48

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

3 of 7  
du  
Bh

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002759  
Core Version: 8105.48  
03/05/2015 13:15:19

---

Standard Lot#: 11514100A2  
Last Changed By: S. CIONE #50616

---

QAS: A. AUTRET #53546  
TUCSON PD

---

Operator: A. AUTRET #53546  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	13:16:21
Diagnostic Test	Pass	13:16:56
Air Blank	0.000	13:17:24
0.100 Cal Check	RFI*	13:17:40
Air Blank	IPS**	13:17:49

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

4 of 7  
db  
B4

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002759  
Core Version: 8105.48  
03/05/2015 13:18:50

---

Standard Lot#: 11514100A2  
Last Changed By: S. CIONE #50616

---

OAS: A. AUTRET #53546  
TUCSON PD

---

Operator: A. AUTRET #53546  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	13:19:51
Diagnostic Test	Pass	13:20:26
Air Blank	0.000	13:20:53
0.100 Cal Check	0.100	13:21:14
Air Blank	0.000	13:21:43
Subject Test	0.080*	13:24:51
Air Blank	0.000	13:25:23
Air Blank	0.000	13:25:50
Subject Test	XXX**	13:26:08
Air Blank	0.000	13:26:42

\*Deficient Sample  
\*\*Invalid Sample  
Begin new deprivation period

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

5 of 7  
MB  
Bh

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002759  
Core Version: 8105.48  
03/05/2015 13:28:04

-----  
Standard Lot#: 11514100A2  
Last Changed By: S. CIONE #50616  
-----

QAS: A. AUTRET #53546  
TUCSON PD  
-----

Operator: A. AUTRET #53546  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	13:28:36
0.100 Cal Check	0.098	13:28:57
Air Blank	0.000	13:29:26

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

6 of 7  
TS4  
3.11.2015

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Amy Autret 53546 AGENCY Tucson Police Department  
DATE 3/5/15 TIME 1306  
Intoxilyzer Serial # 80-002759 Location ODW

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE 

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002759  
Core Version: 8105.48  
03/05/2015

13:06:06

-----  
Standard Lot#: 11514100A2  
Last Changed By: S. CIONE #50616  
-----

QAS: A. AUTRET #53546  
TUCSON PD  
-----

Operator: A. AUTRET #53546  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	13:07:01
0.100 Cal Check	0.099	13:07:22
Air Blank	0.000	13:07:51

Preventative Maintenance  
Performed -- Timer Reset

7 of 7  
ak  
TS4