



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TS4
8-25-2015

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002759 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 8/20/2015 Cal 1 0.096 AC Cal 2 0.095 AC

Between the dates of 7/24/2015 and 8/20/2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on .
The instrument was returned to service on .

Amy L. Autret

Amy L. Autret


8/24/2015

ADPS QAS Permit # 33718

SUBSCRIBED AND SWORN to before me this 24th, day of August, 2015.

Elma Escarcega

Notary Public

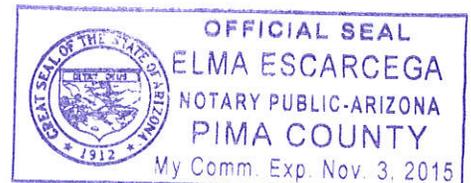

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD QUALITY ASSURANCE PROCEDURE

1087
JLB
TJH
8-29-2015

QA SPECIALIST Amy Autret 53546 AGENCY Tucson Police Department
DATE 8/20/2015 TIME 1502
INTOXILYZER SERIAL # 80-002759 LOCATION ODW

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.


OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 095 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE

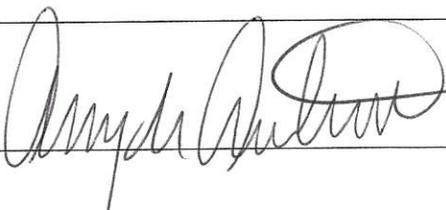

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002759
Core Version: 8105.48
08/20/2015

15:03:32

2017
MA
TS4

Standard Lot#: 11514100A2
Last Changed By: S. CIONE #50616

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 15:04:20 |
| Diagnostic Test | Pass | 15:04:55 |
| Air Blank | 0.000 | 15:05:23 |
| 0.100 Cal Check | 0.096 | 15:05:44 |
| Air Blank | 0.000 | 15:06:13 |
| Subject Test | 0.000 | 15:06:28 |
| Air Blank | 0.000 | 15:06:55 |
| Wait | ABT* | |
| Air Blank | 0.000 | 15:07:29 |

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002759
Core Version: 8105.48
08/20/2015

15:07:48

3 of 7
dd
PH

Standard Lot#: 11514100A2
Last Changed By: S. CIONE #50616

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 15:08:42 |
| Diagnostic Test | Pass | 15:09:17 |
| Air Blank | 0.000 | 15:09:44 |
| 0.100 Cal Check | RFI* | 15:10:00 |
| Air Blank | IPS** | 15:10:09 |

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002759
Core Version: 8105.48
08/20/2015

15:10:18

4 of 7
dlb
TS4

Standard Lot#: 11514100A2
Last Changed By: S. CIONE #50616

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 15:11:10 |
| Diagnostic Test | Pass | 15:11:45 |
| Air Blank | 0.000 | 15:12:12 |
| 0.100 Cal Check | 0.095 | 15:12:33 |
| Air Blank | 0.000 | 15:13:03 |
| Subject Test | 0.296* | 15:16:10 |
| Air Blank | 0.000 | 15:16:44 |
| Air Blank | 0.000 | 15:17:11 |
| Subject Test | XXX** | 15:17:33 |
| Air Blank | 0.000 | 15:18:08 |

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002759
Core Version: 8105.48
08/20/2015 15:18:23

5067
JLB

Standard Lot#: 11514100A2
Last Changed By: S. CIONE #50616

TS4

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 15:18:51 |
| 0.100 Cal Check | 0.095 | 15:19:12 |
| Air Blank | 0.000 | 15:19:42 |

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6 of 7
dlb

TS4
8.25.2015

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Amel Autret 53546 AGENCY Tucson Police Department
DATE 8/20/2015 TIME 1502
Intoxilyzer Serial # 80-002759 Location ODW

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 096 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE 

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002759
Core Version: 8105.48
08/20/2015 15:02:02

7067
Ad
TS4

Standard Lot#: 11514100A2
Last Changed By: S. CIONE #50616

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 15:02:31 |
| 0.100 Cal Check | 0.096 | 15:02:51 |
| Air Blank | 0.000 | 15:03:21 |

Preventative Maintenance
Performed -- Timer Reset