



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 1306 W. Miracle Mile, Tucson, AZ 85705

Tsh
 11-6-2015

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002759 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 10/14/2015 Cal 1 0.096 AC Cal 2 0.095 AC

Between the dates of 9/17/2015 and 10/14/2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on _____.
 The instrument was returned to service on _____.

William L Loomis

William L. Loomis


11/5/2015

ADPS QAS Permit # 33720

SUBSCRIBED AND SWORN to before me this 5th day of November, 2015.

A. M. Hernandez

Notary Public

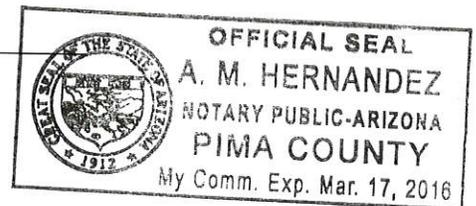

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

1 of 7
WLL
TJG
11-6-2015

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST W. LOOMIS 53545 AGENCY Tucson Police Department

DATE 10/14/15 TIME 17:31

INTOXILYZER SERIAL # 80-002759 LOCATION CRIME LAB

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 095 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE William Loomis

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002759
Core Version: 8105.48
10/14/2015

17:34:00

2 of 7
WLL
TS4

Standard Lot#: 11514100A2
Last Changed By: S. CIONE #50616

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1987
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:34:56
Diagnostic Test	Pass	17:35:31
Air Blank	0.000	17:35:58
0.100 Cal Check	0.096	17:36:19
Air Blank	0.000	17:36:49
Subject Test	0.000	17:37:07
Air Blank	0.000	17:37:34
Wait	ABT*	
Air Blank	0.000	17:38:00

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002759
Core Version: 8105.48
10/14/2015

17:38:32

Standard Lot#: 11514100A2
Last Changed By: S. CIONE #50616

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1987
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:39:21
Diagnostic Test	Pass	17:39:57
Air Blank	0.000	17:40:24
0.100 Cal Check	RFI*	17:40:38
Air Blank	IPS**	17:40:46

*RFI Detect

**Improper Sample

3 of 7 will
TS4

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002759
Core Version: 8105.48
10/14/2015

17:41:07

Standard Lot#: 11514100A2
Last Changed By: S. CIGNE #50616

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1987
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:42:02
Diagnostic Test	Pass	17:42:37
Air Blank	0.000	17:43:04
0.100 Cal Check	0.097	17:43:25
Air Blank	0.000	17:43:54
Subject Test	0.069*	17:47:02
Air Blank	0.000	17:47:34
Air Blank	0.000	17:48:01
Subject Test	XXX**	17:48:21
Air Blank	0.000	17:48:54

*Deficient Sample

**Invalid Sample

Begin new deprivation period

4 of 7
wll
B4

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002759
Core Version: 8105.48
10/14/2015 17:49:44

5 of 7
WLL
TS4

Standard Lot#: 11514100A2
Last Changed By: S. CIONE #50616

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	17:50:14
0.100 Cal Check	0.095	17:50:35
Air Blank	0.000	17:51:04

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7
TS4
11-6-2015

QA Specialist W. LOOMIS 53515 AGENCY Tucson Police Department
DATE 10/14/15 TIME 17:31
Intoxilyzer Serial # 80-002759 Location CRIME LAB

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C.
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 096 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Will Loomis

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002759
Core Version: 8105.48
10/14/2015 17:31:40

7 of 7
WLL
TS4

Standard Lot#: 11514100A2
Last Changed By: S. CIONE #50616

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	17:32:10
0.100 Cal Check	0.096	17:32:31
Air Blank	0.000	17:33:01

Preventative Maintenance
Performed -- Timer Reset