



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002759 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4  
2.11.2016

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 2/4/2016 Cal 1 0.096 AC Cal 2 0.096 AC

Between the dates of 1/7/2016 and 2/4/2016  
the instrument was operating properly and accurately.

The instrument was taken out of service for repair on \_\_\_\_\_.  
The instrument was returned to service on \_\_\_\_\_.

William L Loomis

William L. Loomis

2/10/2016

ADPS QAS Permit # 33720

SUBSCRIBED AND SWORN to before me this 10<sup>th</sup>, day of  
February, 2016.

Elma Escarcega

Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

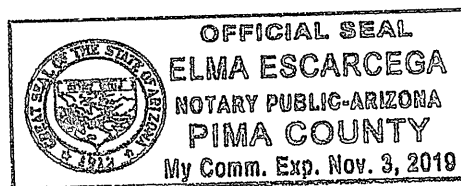


EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7  
WLL  
Bh 2.11.2016

QA SPECIALIST W. LOOMIS 53545 AGENCY Tucson Police Department

DATE 2/4/2016 TIME 15:49

INTOXILYZER SERIAL # 80-002759 LOCATION CRIME LAB

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.  
(✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.  
(✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed  
(✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.  
(✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 096 AC

Instrument is operating properly and accurately. YES ✓ NO     

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE

William Loomis

INTOXILYZER 8000

Location: TUCSON PD

Serial Number: 80-002759

Core Version: 8105.48

02/04/2016

15:51:06

Standard Lot#: 11514100A2

Last Changed By: S. CIONE

#50616

QAS: W. LOOMIS

#53545

TUCSON PD

Operator: W. LOOMIS

#53545

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1987

Sex: M

Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:51:57
Diagnostic Test	Pass	15:52:32
Air Blank	0.000	15:52:59
0.100 Cal Check	0.097	15:53:20
Air Blank	0.000	15:53:50
Subject Test	0.000	15:54:37
Air Blank	0.000	15:55:04
Wait	ABT*	
Air Blank	0.000	15:55:47

\*Sequence Aborted

2 of 7  
WU  
B4

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002759  
Core Version: 8105.48  
02/04/2016

15:59:17

Standard Lot#: 11514100A2  
Last Changed By: S. CIONE #50616

QAS: W. LOOMIS #53545  
TUCSON PD

Operator: W. LOOMIS #53545  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1987  
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:00:29
Diagnostic Test	Pass	16:01:04
Air Blank	0.000	16:01:32
0.100 Cal Check	RFI*	16:01:45
Air Blank	IPS**	16:01:54

\*RFI Detect

\*\*Improper Sample

3 of 7  
WLL  
PL

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000

Location: TUCSON PD

Serial Number: 80-002759

Core Version: 8105.48

02/04/2016

16:02:01

Standard Lot#: 11514100A2

Last Changed By: S. CIONE

#50616

QAS: W. LOOMIS

#53545

TUCSON PD

Operator: W. LOOMIS

#53545

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1987

Sex: M

Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:02:52
Diagnostic Test	Pass	16:03:27
Air Blank	0.000	16:03:55
0.100 Cal Check	0.097	16:04:16
Air Blank	0.000	16:04:45
Subject Test	0.252*	16:07:53
Air Blank	0.000	16:08:25
Air Blank	0.000	16:08:53
Subject Test	XXX**	16:09:11
Air Blank	0.000	16:09:45

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

4 of 7  
WLL  
B4

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000

Location: TUCSON PD

Serial Number: 80-002759

Core Version: 8105.48

02/04/2016

16:13:11

-----  
Standard Lot#: 11514100A2

Last Changed By: S. CIONE #50616  
-----

QAS: W. LOOMIS #53545

TUCSON PD  
-----

Operator: W. LOOMIS #53545

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	16:13:41
0.100 Cal Check	0.096	16:14:01
Air Blank	0.000	16:14:31

Preventative Maintenance  
Performed -- Timer Reset

5 of 7  
WLL  
TS4

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7  
WLL  
By 2.11.2016

QA Specialist W. LOOMIS 53545 AGENCY Tucson Police Department  
DATE 2/4/2016 TIME 15:49  
Intoxilyzer Serial # 80-002759 Location CRIME LAB

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 096 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE William Loomis

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON PD

Serial Number: 80-002759

Core Version: 8105.48

02/04/2016

15:49:34

-----  
Standard Lot#: 11514100A2

Last Changed By: S. CIONE #50616  
-----

QAS: W. LOOMIS #53545

TUCSON PD  
-----

Operator: W. LOOMIS #53545

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	15:50:03
0.100 Cal Check	0.096	15:50:24
Air Blank	0.000	15:50:54

Preventative Maintenance  
Performed -- Timer Reset

7 of 7  
WLL  
B4