



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

TS4  
8.24.2016

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002759 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 8/23/2016 Cal 1 0.097 AC Cal 2 0.096 AC

Between the dates of 7/21/2016 and 8/23/2016 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on \_\_\_\_\_.  
The instrument was returned to service on \_\_\_\_\_.

William L Loomis

William L. Loomis


8/24/2016

ADPS QAS Permit # 33720

SUBSCRIBED AND SWORN to before me this 24<sup>th</sup> day of August, 2016.

Elma Escarcega

Notary Public

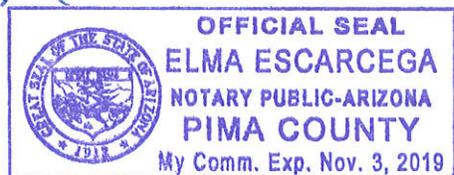

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7  
WLV  
TS4 8.24.2016

QA SPECIALIST W. LOOMIS 53545 AGENCY Tucson Police Department  
DATE 8/23/2016 TIME 15:34  
INTOXILYZER SERIAL # 80-002759 LOCATION CRIME LAB

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 096 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS 731 days  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Will Loomis

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002759  
Core Version: 8105.48  
08/23/2016

15:37:53

Standard Lot#: 06616100A1  
Last Changed By: W. LOOMIS #53545

QAS: W. LOOMIS #53545  
TUCSON PD

Operator: W. LOOMIS #53545  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1987  
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:38:50
Diagnostic Test	Pass	15:39:25
Air Blank	0.000	15:39:52
0.100 Cal Check	0.097	15:40:13
Air Blank	0.000	15:40:42
Subject Test	0.000	15:40:58
Air Blank	0.000	15:41:26
Wait	ABT*	
Air Blank	0.000	15:42:01

\*Sequence Aborted

2 of 7  
WLL  
TS4

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002759  
Core Version: 8105.48  
08/23/2016 15:42:06

Standard Lot#: 06616100A1  
Last Changed By: W. LOOMIS #53545

QAS: W. LOOMIS #53545  
TUCSON PD

Operator: W. LOOMIS #53545  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1987  
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:42:55
Diagnostic Test	Pass	15:43:30
Air Blank	0.000	15:43:57
0.100 Cal Check	RFI*	15:44:11
Air Blank	IPS**	15:44:19

\*RFI Detect  
\*\*Improper Sample

3 of 7  
WLL  
TS4

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002759  
Core Version: 8105.48  
08/23/2016

15:45:56

Standard Lot#: 06616100A1  
Last Changed By: W. LOOMIS #53545

QAS: W. LOOMIS #53545  
TUCSON PD

Operator: W. LOOMIS #53545  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1987  
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:46:43
Diagnostic Test	Pass	15:47:18
Air Blank	0.000	15:47:46
0.100 Cal Check	0.097	15:48:07
Air Blank	0.000	15:48:36
Subject Test	0.103*	15:51:43
Air Blank	0.000	15:52:16
Air Blank	0.000	15:52:43
Subject Test	XXX**	15:53:02
Air Blank	0.000	15:53:35

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

4 of 7  
WLL  
TS4

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002759  
Core Version: 8105.48  
08/23/2016 15:54:08

-----  
Standard Lot#: 06616100A1  
Last Changed By: W. LOOMIS #53545  
-----

QAS: W. LOOMIS #53545  
TUCSON PD  
-----

Operator: W. LOOMIS #53545  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	15:54:38
0.100 Cal Check	0.096	15:54:59
Air Blank	0.000	15:55:29

Preventative Maintenance  
Performed -- Timer Reset

5 of 7  
b WLL  
TSh

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7  
WLL  
TS 8.24.2016

QA Specialist W. LOOMIS 53545 AGENCY Tucson Police Department  
DATE 8/23/2016 TIME 15:34  
Intoxilyzer Serial # 80-002759 Location CRIME LAB

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads ~~"PUSH BUTTON TO START"~~ <sup>CONTACT GAS- 31 DAY PREV MAINT</sup>
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 097 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Willie Loomis

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002759  
Core Version: 8105.48  
08/23/2016 15:34:28

-----  
Standard Lot#: 06616100A1  
Last Changed By: W. LOOMIS #53545  
-----

QAS: W. LOOMIS #53545  
TUCSON PD  
-----

Operator: W. LOOMIS #53545  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	15:34:59
0.100 Cal Check	0.097	15:35:20
Air Blank	0.000	15:35:49

Preventative Maintenance  
Performed -- Timer Reset

7 of 7  
WLL  
TS4