



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

cu
11.23.07

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002760 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 10.23.2007 Cal 1 0.103 AC Cal 2 0.101 AC

The instrument was placed into service on 10.23.2007.


Terry S. Gallegos

Terry S. Gallegos

11.21.2007

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 21, day of November, 2007.

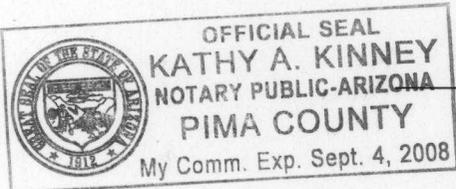

Kathy A. Kinney

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

CU
11-23-07

CERTIFICATE OF CALIBRATION

INTOXILYZER® 8000

This is to certify that the calibration of INTOXILYZER serial number 80-002760, manufactured by CMI, subsidiary of MPD, Inc. of Owensboro, Kentucky, was tested and found to be in compliance with the National Safety Administration Standard for Devices to Measure Breath Alcohol (F.R., No. 179 48705-48710 Sept. 17, 1993) and with International Recommendation OIML R126 Edition 1998 (E). Calibration solutions are traceable to NIST (NBS) standard material 1828.

Date 10-5-07

Signed Dale Howard
Technician

CMI INC.

316 EAST 9TH STREET
OWENSBORO, KY. 42303



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

CM
11.23.07

INTOXILYZER SERVICE RECORD

10/1
T84

Intoxilyzer Serial Number: 80-002760

Initial Service Date: 10-23-2007

CMI, Inc. Certificate of Calibration Date: 10.5.2007

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

cu
11-23-07

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

1 of 7
TS4

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 10.23.2007 TIME 1750

INTOXILYZER SERIAL # 80-002760 LOCATION Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS


- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed ✓
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 101 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

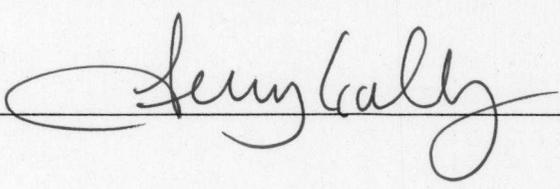SIGNATURE 

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002760
10/23/2007

17:52:31

CM
11-23-07
2 of 7
TB4

Standard Lot#: 725401I
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:53:21
Diagnostic Test	Pass	17:53:49
Air Blank	0.000	17:54:16
0.100 Cal Check	0.102	17:54:35
Air Blank	0.000	17:55:05
Subject Test	0.000	17:56:02
Air Blank	0.000	17:56:30
Five Minute Wait	ABT*	
Air Blank	0.000	17:57:09

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002760
10/23/2007

17:57:40

CU
11-23-07

3 of 7
TS4

Standard Lot#: 725401I
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:58:29
Diagnostic Test	Pass	17:58:56
Air Blank	0.000	17:59:23
0.100 Cal Check	0.103	17:59:42
Air Blank	0.000	18:00:12
Subject Test	RFI*	18:00:52
Air Blank	IPS**	18:01:00

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1527, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002760
10/23/2007 18:01:08

CU
11-23-07
4 of 7
TSh

Standard Lot#: 725401I
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	18:01:56
Diagnostic Test	Pass	18:02:24
Air Blank	0.000	18:02:51
0.100 Cal Check	0.102	18:03:10
Air Blank	0.000	18:03:39
Subject Test	0.101*	18:06:48
Air Blank	0.000	18:07:20
Air Blank	0.000	18:07:48
Subject Test	XXX**	18:08:14
Air Blank	0.000	18:08:48

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1527, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002760
10/23/2007 18:09:06

en

Soft
TSh

Standard Lot#: 725401I
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	18:09:33
Cal Check	0.101	18:09:52
Air Blank	0.000	18:10:21

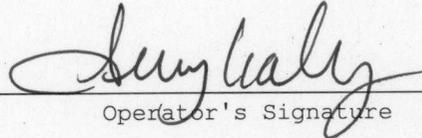

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

CM 23-07
11

ARIZONA DEPARTMENT OF PUBLIC SAFETY


STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

6077
TSH

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28947 AGENCY Tucson Police Department
DATE 10.23.2007 TIME 1750
Intoxilyzer Serial # 80-002760 Location Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 103 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

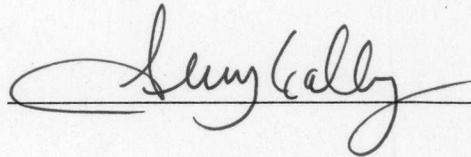SIGNATURE 

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002760
10/23/2007 17:50:34

cm

7 of 7
B4

Standard Lot#: 725401I
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	17:51:01
Cal Check	0.103	17:51:20
Air Blank	0.000	17:51:50



Operator's Signature