



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

704
3.17.08

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002760 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3-7-08 Cal 1 0.101 AC Cal 2 0.101 AC

Between the dates of 2-8-08 and 3-7-08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.


Clare Conley

Clare L. Conley

3-14-08

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 14, day of March, 2008.

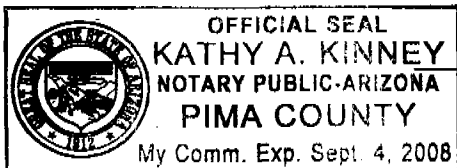

Kathy A Kinney
Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Clare Conley AGENCY Tucson Police Department
DATE 3-7-08 TIME 0644
INTOXILYZER SERIAL # 80-002760 LOCATION ODE

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0.000 AC.
(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(✓) 4. Standard Calibration Check standard 0.100 AC
Result: 0.101 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS None.

SIGNATURE Clare Conley

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002760
03/07/2008

06:46:46

2 of 7
Cn
154

Standard Lot#: 725401I
Last Changed By: T. GALLEGOS #28847
QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	06:47:37
Diagnostic Test	Pass	06:48:04
Air Blank	0.000	06:48:32
0.100 Cal Check	0.101	06:48:50
Air Blank	0.000	06:49:20
Subject Test	0.000	06:49:35
Air Blank	ABT*	06:49:41
Air Blank	0.000	06:50:10

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002760
03/07/2008

06:50:35

3:17
cu
TS4

Standard Lot#: 7254011
Last Changed By: T. GALLEGOS #28847

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TESI, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	06:51:17
Air Blank	IPS**	06:51:26

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002760
03/07/2008

06:51:50

Standard Lot#: 725401I

Last Changed By: T. GALLEGOS #28847

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST,SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	06:52:47
Diagnostic Test	Pass	06:53:14
Air Blank	0.000	06:53:41
0.100 Cal Check	0.101	06:54:00
Air Blank	0.000	06:54:30
Subject Test	0.107*	06:57:38
Air Blank	0.000	06:58:10
Air Blank	0.000	06:58:38
Subject Test	XXX**	06:58:59
Air Blank	0.000	06:59:33

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002760
03/07/2008 07:00:04

55f7
An
TS4

Standard Lot#: 7254011
Last Changed By: T. GALLEGOS #28847

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	07:00:31
Cal Check	0.101	07:00:50
Air Blank	0.000	07:01:20

Clare Conley
Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 3-7-08 TIME 0644
Intoxilyzer Serial # 80-002760 Location DDE

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002760
03/07/2008

06:44:27

Standard Lot#: 7254011
Last Changed By: T. GALLEGOS #28847

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	06:44:55
Cal Check	0.101	06:45:14
Air Blank	0.000	06:45:43

Clare Conley

Operator's Signature