



CERTIFICATION by CUSTODIAN of RECORDS  
Tucson Police Department Crime Laboratory

By  
6-30-08

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002760 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6-27-08 Cal 1 0.099 AC Cal 2 0.098 AC

Between the dates of 5-30-08 and 6-27-08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A  
The instrument was returned to service on N/A


Clare Conley

Clare L. Conley

6-27-08

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 27, day of June, 2008.

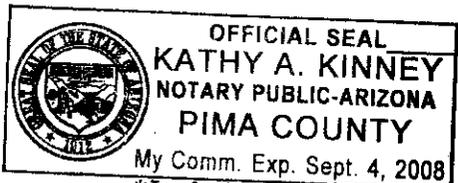

Kathy A Kinney

Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1 of 7  
cu

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

By  
6-30-08

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 6-27-08 TIME 0642

INTOXILYZER SERIAL # 80-002760 LOCATION ODE

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 098 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS None  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Clare Conley

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002760  
06/27/2008

06:44:47

2007  
ew  
TS4

Standard Lot#: 725401I  
Last Changed By: T. GALLEGOS #28847

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	06:45:40
Diagnostic Test	Pass	06:46:07
Air Blank	0.000	06:46:35
0.100 Cal Check	0.100	06:46:54
Air Blank	0.000	06:47:23
Subject Test	0.000	06:47:46
Air Blank	ABT*	06:47:53
Air Blank	0.000	06:48:21

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002760  
06/27/2008

06:48:26

Standard Lot#: 725401I  
Last Changed By: T. GALLEGOS #28847

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	06:48:59
Air Blank	IPS**	06:49:07

\*RFI Detect

\*\*Improper Sample

387  
04  
104

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002760  
06/27/2008

06:49:20

Standard Lot#: 725401I  
Last Changed By: T. GALLEGOS #28847

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	06:50:10
Diagnostic Test	Pass	06:50:37
Air Blank	0.000	06:51:05
0.100 Cal Check	0.099	06:51:23
Air Blank	0.000	06:51:53
Subject Test	0.185*	06:55:01
Air Blank	0.000	06:55:33
Air Blank	0.000	06:56:01
Subject Test	XXX**	06:56:29
Air Blank	0.000	06:57:03

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

457  
Cu  
BY

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002760  
06/27/2008 06:57:31

ESB7  
CW  
TS4

-----  
Standard Lot#: 725401I  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: C. CONLEY #34517  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	06:57:58
Cal Check	0.098	06:58:17
Air Blank	0.000	06:58:46

*Clare Conley*

\_\_\_\_\_  
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6-27-08  
154  
6-30-08  
CW

QA Specialist Clare Conley AGENCY Tucson Police Department  
DATE 6-27-08 TIME 0642  
Intoxilyzer Serial # 80-002760 Location ODE

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0.100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach  $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0.099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002760  
06/27/2008 06:42:55

-----  
Standard Lot#: 725401I  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: C. CONLEY #34517  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	06:43:22
Cal Check	0.099	06:43:41
Air Blank	0.000	06:44:11

*Clare Conley*

Operator's Signature

787  
754  
cu