



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002760 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4
11-4-2008

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 10-17-08 Cal 1 0.099 AC Cal 2 0.097 AC

Between the dates of 9-19-08 and 10-17-08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/APP.

Clare Conley

Clare L. Conley


11-4-08

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 4th day of November, 2008.

Ant Hernandez

Notary Public

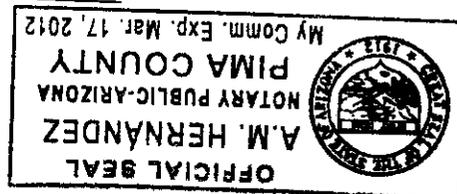

*Includes documents specified in A.A.C. R13-10-104.A.5

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1067
Cu
TS4
11-4-2008

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 10-17-08 TIME 0733

INTOXILYZER SERIAL # 80-002760 LOCATION DDE

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.

(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.

(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed

(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.

(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 097 AC

Instrument is operating properly and accurately. YES NO

COMMENTS None

SIGNATURE Clare Conley

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002760
10/17/2008

07:33:24

2 of 7
Cu
104

Standard Lot#: 725401I
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	07:34:08
Air Blank	IPS**	07:34:24

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002760
10/17/2008

07:36:55

387
en
154

Standard Lot#: 725401I
Last Changed By: T. GALLEGOS #28847

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:37:48
Diagnostic Test	Pass	07:38:15
Air Blank	0.000	07:38:43
0.100 Cal Check	0.098	07:39:02
Air Blank	0.000	07:39:31
Subject Test	0.000	07:39:47
Air Blank	ABT*	07:39:53
Air Blank	0.000	07:40:22

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002760
10/17/2008

07:40:26

4367
lu
TS4

Standard Lot#: 725401I
Last Changed By: T. GALLEGOS #28847
OAS: C. CONLEY #34517
TUCSON PD
Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SOAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:41:12
Diagnostic Test	Pass	07:41:39
Air Blank	0.000	07:42:07
0.100 Cal Check	0.099	07:42:25
Air Blank	0.000	07:42:55
Subject Test	0.466*	07:46:04
Air Blank	0.000	07:46:37
Air Blank	0.000	07:47:05
Subject Test	XXX**	07:47:27
Air Blank	0.000	07:48:02

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002760
10/17/2008 07:48:19

5687
lu

Standard Lot#: 725401I
Last Changed By: T. GALLEGOS #28847

154

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	07:48:46
Cal Check	0.097	07:49:05
Air Blank	0.000	07:49:35

Clare Conley
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD CALIBRATION CHECK PROCEDURE

6587
Cu
TS4
11.4.2008

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 10-17-08 TIME 0733
Intoxilyzer Serial # 80-002760 Location ODE

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 90-002760
10/17/2008

07:34:55

7587
CW
154

Standard Lot#: 725401I
Last Changed By: T. GALLEGOS #28847

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	07:35:22
Cal Check	0.099	07:35:41
Air Blank	0.000	07:36:10

Clare Conley
Operator's Signature