



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

CM
11-21-08

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002760 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 11/14/08 Cal 1 0.099 AC Cal 2 0.098 AC

Between the dates of 10/17/08 and 11/14/08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on _____.

The instrument was returned to service on _____.

Sarah Cione

Sarah L. Cione


11/20/08

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 20th day of November, 2008.

A. M. Hernandez

Notary Public

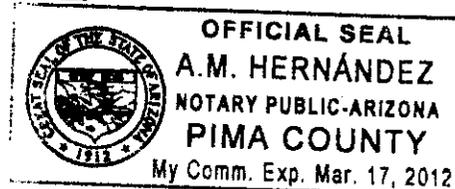

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 10_{5c}

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

Ch
11-21-08

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 11/14/08 TIME 10:20

INTOXILYZER SERIAL # 80-002760 LOCATION Team 4

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 098 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

2 of 10 SC
W

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002760
11/14/2008

10:21:38

Standard Lot#: 725401I
Last Changed By: T.- GALLEGOS #28847
QAS: S. CIONE #50616
TUCSON PD
Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:22:25
Diagnostic Test	Pass	10:22:52
Air Blank	0.000	10:23:20
0.100 Cal Check	0.099	10:23:38
Air Blank	0.000	10:24:08
Subject Test	0.000	10:24:27
Air Blank	0.000	10:24:54
Five Minute Wait	ABT*	
Air Blank	0.000	10:25:37

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

3 of 10 SC
cu

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002760
11/14/2008

10:25:50

Standard Lot#: 725401I
Last Changed By: T.-GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:26:54
Diagnostic Test	Pass	10:27:21
Air Blank	0.000	10:27:49
0.100 Cal Check	0.099	10:28:08
Air Blank	0.000	10:28:37
Subject Test	0.100*	10:31:45
Air Blank	0.000	10:32:18
Air Blank	0.000	10:32:45
Subject Test	XXX**	10:33:09
Air Blank	0.000	10:33:42

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

4 of 10 SC

CM

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002760
11/14/2008

10:33:48

Standard Lot#: 725401I
Last Changed By: T.-GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	10:34:14
Air Blank	IPS**	10:34:24

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

5 of 10 SC

cu

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002760
11/14/2008 10:34:45

Standard Lot#: 725401I
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	AMB*	10:35:13
Air Blank	0.000	10:35:41

*Ambient Fail

Sarah Cione
Operator's Signature

performed after
improper sample,
of which there was
still breath spray
present in my mouth

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002760
11/14/2008

10:35:48

b of 10 sc

cu

Standard Lot#: 725401I
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:36:15
Cal Check	0.098	10:36:34
Air Blank	0.000	10:37:03

Sarah Cione

Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

7 of 10 SC
CW
11-21-08

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 11/14/08 TIME 10:20
Intoxilyzer Serial # 80-002760 Location Team 4

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002760
11/14/2008

10:20:13

8 of 10 SC

cu

Standard Lot#: 725401I
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:20:40
Cal Check	0.099	10:20:59
Air Blank	0.000	10:21:29

Sarah Cione

Operator's Signature



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

9 of 10_{sc}

Ch
11-21-08

INTOXILYZER SERVICE RECORD
FOR
CALIBRATION GAS TANK CHANGE

Date: 11/14/08
 QAS: Sarah Cione
 Intoxilyzer Serial Number: 80-002760
 New Lot Number: 808101G
 New Expiration Date: 3/31/2010

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change." Select "D" for dry gas and enter the following information:
 Enter Target Value: 0.100
 Enter Control Lot #: From the tank label
 Enter Exp Date: From the tank label (mm/dd/yyyy)
 Edit info Y/N? Y or N
 Update Std Info? Y ENTER
 "Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing." Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 102 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove and sign printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002760
11/14/2008

10:39:18

10 of 10 SC

cu

Standard Lot#: 808101G
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:39:45
Cal Check	0.102	10:40:04
Air Blank	0.000	10:40:34

Sarah Cione

Operator's Signature