



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

TS4  
3-17-09

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002760 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3/6/09 Cal 1 0.101 AC Cal 2 0.101 AC

Between the dates of 2/6/09 and 3/6/09 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -  
The instrument was returned to service on -

Sarah Cione

Sarah L. Cione


3/17/09

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 17, day of March, 2009.

Art Hernandez

Notary Public

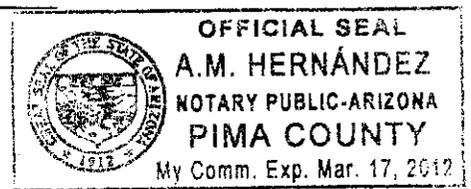

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

1 of 7 SC

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD QUALITY ASSURANCE PROCEDURE

TS4  
3-17-09

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department  
DATE 3/6/09 TIME 12:08  
INTOXILYZER SERIAL # 80-002760 LOCATION Team 4

( / ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 101 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2 of 75c

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002760  
03/06/2009

12:10:04

154

Standard Lot#: 808101G  
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP,TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:10:50
Diagnostic Test	Pass	12:11:17
Air Blank	0.000	12:11:45
0.100 Cal Check	0.101	12:12:03
Air Blank	0.000	12:12:33
Subject Test	0.000	12:12:49
Air Blank	0.000	12:13:16
Five Minute Wait	ABT*	
Air Blank	0.000	12:13:52

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

3 of 7 SC

154

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002760  
03/06/2009

12:13:56

Standard Lot#: 808101G  
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	12:14:28
Air Blank	IPS**	12:14:36

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

4 of 75C  
154

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002760  
03/06/2009 12:14:42  
Standard Lot#: 808101G  
Last Changed By: S. CIONE #50616  
QAS: S. CIONE #50616  
TUCSON PD  
Operator: S. CIONE #50616  
TUCSON PD  
Subject: SQAP, TEST  
DOB: 11/11/1111  
Sex: F Weight: 111  
15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:15:26
Diagnostic Test	Pass	12:15:53
Air Blank	0.000	12:16:20
0.100 Cal Check	0.100	12:16:39
Air Blank	0.000	12:17:09
Subject Test	0.097*	12:20:17
Air Blank	0.000	12:20:50
Air Blank	0.000	12:21:17
Subject Test	XXX**	12:21:37
Air Blank	0.000	12:22:11

\*Deficient Sample  
\*\*Invalid Sample  
Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 86-062760  
03/05/2009

12:22:27

5 of 7 SC

Standard Lot#: 0081010  
Last Changed By: S. CIONE #50616

TS4

OAS: S. CIONE #50616  
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	12:22:54
Cal Check	0.101	12:23:13
Air Blank	0.000	12:23:43

*Sarah Cione*

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

6 of 7 SC

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

TS4  
3-17-09

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 3/6/09 TIME 12:08  
Intoxilyzer Serial # 80-002760 Location Team 4

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002760  
03/06/2009 12:08:59

7 of 7 SC  
TSY

-----  
Standard Lot#: 898101G  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	12:09:07
Cal Check	0.101	12:09:25
Air Blank	0.000	12:09:55

*Sarah Cione*  
-----  
Operator's Signature