



CERTIFICATION by CUSTODIAN of RECORDS  
Tucson Police Department Crime Laboratory

9.2.09

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002760 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 8.20.2009 Cal 1 0.098 AC Cal 2 0.097 AC

Between the dates of 7.23.2009 and 8.20.2009 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on                     .  
The instrument was returned to service on                     .

Terry S. Gallegos

Terry S. Gallegos


9.2.2009

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 2nd, day of September, 2009.

A.M. Hernandez

Notary Public

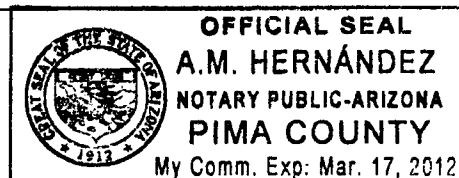

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

CU  
9.2-09

1 of 7  
TS6

QA SPECIALIST Terry Gallegos 20047 AGENCY Tucson Police Department  
DATE 8.20.2009 TIME 1144  
INTOXILYZER SERIAL # 80-002760 LOCATION DUI HQ

( ☒ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ☒ ) 1. Clock time check.  
( ☒ ) 2. Date check.

OPERATIONAL TESTS

- ( ☒ ) 1. Alcohol-free subject test result 0. 000 AC.  
( ☒ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed ✓✓✓  
( ☒ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.  
( ☒ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 097 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE

Terry Gallegos

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002760  
08/20/2009

11:46:23

Standard Lot#: 730901I  
Last Changed By: S. CIONE #50616

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:47:33
Diagnostic Test	Pass	11:48:01
Air Blank	0.000	11:48:28
0.100 Cal Check	0.098	11:48:47
Air Blank	0.000	11:49:16
Subject Test	0.000	11:50:02
Air Blank	ABT*	11:50:14
Air Blank	0.000	11:50:42

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

*Pursuant to section 28-1527, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.*

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002760  
08/20/2009

11:50:52

Standard Lot#: 730901I  
Last Changed By: S. CIONE #50616  
QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:51:58
Diagnostic Test	Pass	11:52:26
Air Blank	RFI*	11:52:38
Air Blank	IPS**	11:52:46

\*RFI Detect

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

*Pursuant to section 28-1527, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department*

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002760  
08/20/2009

11:53:12

Standard Lot#: 730901I  
Last Changed By: S. CIONE #50616

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:54:02
Diagnostic Test	Pass	11:54:30
Air Blank	0.000	11:54:57
0.100 Cal Check	0.098	11:55:16
Air Blank	0.000	11:55:45
Subject Test	0.266*	11:58:54
Air Blank	0.000	11:59:26
Air Blank	0.000	11:59:54
Subject Test	XXX**	12:00:35
Air Blank	0.000	12:01:09

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

*Pursuant to section 28-1527, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.*

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002760  
08/20/2009 12:01:28

-----  
Standard Lot#: 730901I  
Last Changed By: S. CIONE #50616  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB  
-----

Test	g/210L	Time
Air Blank	0.000	12:01:56
Cal Check	0.097	12:02:14
Air Blank	0.000	12:02:44

  
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

CM  
a-2-09  
6087  
TJG

QA Specialist Terry Gallegos 20847 AGENCY Tucson Police Department  
DATE 8-20-2009 TIME 1144  
Intoxilyzer Serial # 80-002760 Location DUI HQ

- ( ☒ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach  $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- ( ☒ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ☒ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ☒ ) 4. Air blank completed.
- ( ☒ ) 5. Calibration check completed. Test results 0. 098 AC.
- ( ☒ ) 6. Air blank completed.
- ( ☒ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)


INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002760  
08/20/2009

11:44:24

-----  
Standard Lot#: 730901I  
Last Changed By: S. CIONE #50616  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB  
-----

Test	g/210L	Time
Air Blank	0.000	11:44:52
Cal Check	0.098	11:45:11
Air Blank	0.000	11:45:40

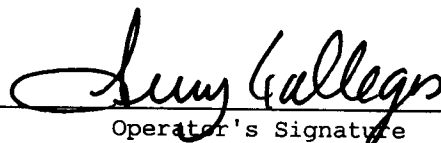  
Operator's Signature