



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

TS4  
5-5-2010

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002760 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 4/29/10 Cal 1 0.101 AC Cal 2 0.100 AC

Between the dates of 4/1/10 and 4/29/10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


4/30/10

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 30<sup>th</sup> day of April, 2010.

A.M. Hernandez

Notary Public

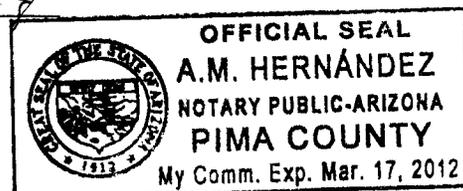

\*Includes documents specified in A.A.C. R13-10-104.A.5

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1 of 7 SC

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

T54  
5.5.2010

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 4/29/10 TIME 858

INTOXILYZER SERIAL # 80-002760 LOCATION DUI HQ - C12

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 100 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE Sarah Cione

2 of 7 sc

B4

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002760  
04/29/2010 09:01:51

Standard Lot#: 826001I  
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test             | g/210L | Time     |
|------------------|--------|----------|
| Air Blank        | 0.000  | 09:02:38 |
| Diagnostic Test  | Pass   | 09:03:05 |
| Air Blank        | 0.000  | 09:03:32 |
| 0.100 Cal Check  | 0.101  | 09:03:51 |
| Air Blank        | 0.000  | 09:04:21 |
| Subject Test     | 0.000  | 09:04:48 |
| Air Blank        | 0.000  | 09:05:16 |
| Five Minute Wait | ABT*   |          |
| Air Blank        | 0.000  | 09:05:57 |

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

3 of 7 sc

BH

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002760  
04/29/2010 09:07:02

Standard Lot#: 826001I  
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 09:07:45 |
| Diagnostic Test | Pass   | 09:08:12 |
| Air Blank       | RFI*   | 09:08:38 |
| Air Blank       | IPS**  | 09:08:52 |

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

*Pursuant to section 28-1527, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.*

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002760  
04/29/2010

09:09:01

4 of 7 sc

124

---

Standard Lot#: 826001I  
Last Changed By: C. CONLEY #34517

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 09:09:59 |
| Diagnostic Test | Pass   | 09:10:27 |
| Air Blank       | 0.000  | 09:10:54 |
| 0.100 Cal Check | 0.100  | 09:11:13 |
| Air Blank       | 0.000  | 09:11:43 |
| Subject Test    | 0.062* | 09:14:51 |
| Air Blank       | 0.000  | 09:15:23 |
| Air Blank       | 0.000  | 09:15:50 |
| Subject Test    | XXX**  | 09:16:44 |
| Air Blank       | 0.000  | 09:17:17 |

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002760  
04/29/2010

09:17:43

5 of 7 sc

-----  
Standard Lot#: 826001I  
Last Changed By: C. CONLEY #34517  
-----

BS4

-----  
QAS: S. CIONE #50616  
TUCSON PD  
-----

| Test      | g/210L | Time     |
|-----------|--------|----------|
| Air Blank | 0.000  | 09:18:11 |
| Cal Check | 0.100  | 09:18:29 |
| Air Blank | 0.000  | 09:18:59 |

*Sarah Cione*

\_\_\_\_\_  
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

6 of 7 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

TS4  
5.5.2010

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 4/29/10 TIME 858  
Intoxilyzer Serial # 80-002760 Location DUI HQ - C12

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002760  
04/29/2010

08:58:47

7 of 7 sc  
BH

-----  
Standard Lot#: 826001I  
Last Changed By: C. CONLEY #34517  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

| Test      | g/210L | Time     |
|-----------|--------|----------|
| Air Blank | 0.000  | 08:59:15 |
| Cal Check | 0.101  | 08:59:33 |
| Air Blank | 0.000  | 09:00:03 |

Sarah Cione  
Operator's Signature