



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

TS4
726-2010

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002760 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 7/22/10 Cal 1 0.101 AC Cal 2 0.100 AC

Between the dates of 6/24/10 and 7/22/10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione

7/22/10


ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 22, day of July, 2010.

A.M. Hernández

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

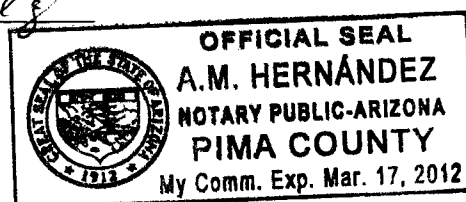

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 7 SC

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000B4
7-26-2010

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police DepartmentDATE 7/22/10 TIME 855INTOXILYZER SERIAL # 80-002760 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS

SIGNATURE Sarah Cione

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002760
07/22/2010

08:58:48

Standard Lot#: 826001I
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:59:51
Diagnostic Test	Pass	09:00:19
Air Blank	0.000	09:00:46
0.100 Cal Check	0.101	09:01:05
Air Blank	0.000	09:01:34
Subject Test	0.000	09:01:56
Air Blank	0.000	09:02:23
Five Minute Wait	ABT*	
Air Blank	0.000	09:03:07

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002760
07/22/2010

09:04:07

3 of 7_{SC}

BH

Standard Lot#: 826001I
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	09:04:40
Air Blank	IPS**	09:04:47

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

4 of 7_{SC}

BY

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002760
07/22/2010

09:05:20

Standard Lot#: 826001I
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:06:27
Diagnostic Test	Pass	09:06:55
Air Blank	0.000	09:07:22
0.100 Cal Check	0.100	09:07:41
Air Blank	0.000	09:08:11
Subject Test	0.041*	09:11:19
Air Blank	0.000	09:11:50
Air Blank	0.000	09:12:18
Subject Test	XXX**	09:12:33
Air Blank	0.000	09:13:06

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002760
07/22/2010

09:13:25

Standard Lot#: 826001I
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:13:53
Cal Check	0.100	09:14:12
Air Blank	0.000	09:14:41

Sarah Cione
Operator's Signature

EXHIBIT G-2

6 of 7_{sc}

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

TS4
7.26.2010

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 7/22/10 TIME 855
Intoxilyzer Serial # 80-002760 Location DUT HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002760
07/22/2010

08:55:31

Standard Lot#: 826001I
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:55:58
Cal Check	0.101	08:56:17
Air Blank	0.000	08:56:46

Sarah Cione
Operator's Signature