



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002760 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4
5.2.2011

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 4/28/11 Cal 1 0.102 AC Cal 2 0.102 AC

Between the dates of 3/31/11 and 4/28/11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -

The instrument was returned to service on -

Sarah Cione

Sarah L. Cione


4/28/11

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 28th day of April, 2011.

A.M. Hernandez

Notary Public

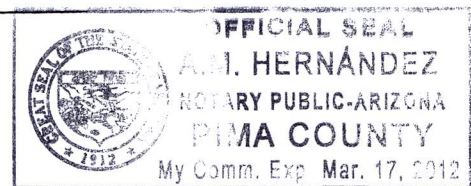

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7 SC

B4
5-2-2011QA SPECIALIST Sarah Cione AGENCY Tucson Police DepartmentDATE 4/28/11 TIME 927INTOXILYZER SERIAL # 80-002760 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 102 AC

Instrument is operating properly and accurately. YES ✓ NO

COMMENTS

SIGNATURE

Sarah Cione

2 of 7_{sc}

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002760
04/28/2011

09:28:53

SM

Standard Lot#: 719919
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:29:39
Diagnostic Test	Pass	09:30:07
Air Blank	0.000	09:30:34
0.100 Cal Check	0.102	09:30:53
Air Blank	0.000	09:31:23
Subject Test	0.000	09:32:41
Air Blank	0.000	09:33:08
Five Minute Wait	ABT*	
Air Blank	0.000	09:33:50

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002760
04/28/2011

09:35:26

3 of 7_{sc}

by

Standard Lot#: 719919
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:36:17
Diagnostic Test	Pass	09:36:45
Air Blank	RFI*	09:37:02
Air Blank	IPS**	09:37:12

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002760
04/28/2011

09:41:55

4 of 7_{sc}
154

Standard Lot#: 719919
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:43:00
Diagnostic Test	Pass	09:43:28
Air Blank	0.000	09:43:55
0.100 Cal Check	0.103	09:44:14
Air Blank	0.000	09:44:44
Subject Test	0.135*	09:47:52
Air Blank	0.000	09:48:25
Air Blank	0.000	09:48:52
Subject Test	XXX**	09:49:08
Air Blank	0.000	09:49:42

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002760
04/28/2011 09:50:22

5 of 7_{sc}

Standard Lot#: 719919
Last Changed By: S. CIONE #50616

TS4

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:50:50
Cal Check	0.102	09:51:08
Air Blank	0.000	09:51:38

Sarah Cione

Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

134
5.2.2011

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 4/28/11 TIME 927
Intoxilyzer Serial # 80-002760 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 102 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

7 of 7
sc

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002760
04/28/2011 09:27:14

Standard Lot#: 719919
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:27:42
Cal Check	0.102	09:28:01
Air Blank	0.000	09:28:30

Sarah Cione
Operator's Signature