



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002760 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4
6-27-2011

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6/23/11 Cal 1 0.101 AC Cal 2 0.101 AC

Between the dates of 5/26/11 and 6/23/11
the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -
The instrument was returned to service on -

Sarah Cione

Sarah L. Cione

6/27/11


ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 27, day of June, 20 11.

A.M. Hernandez

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

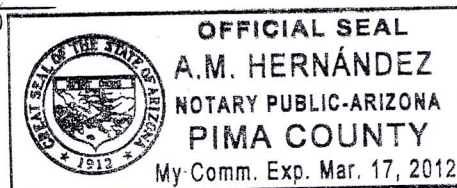

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 9_{sc}

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TSC
6-27-2011

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 6/23/11 TIME 957

INTOXILYZER SERIAL # 80-002760 LOCATION DUI HQ

(☒) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (☒) 1. Clock time check.
(☒) 2. Date check.

OPERATIONAL TESTS

- (☒) 1. Alcohol-free subject test result 0. 000 AC.
(☒) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(☒) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(☒) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 101 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS

SIGNATURE

Sarah Cione

DPS Form Exh G-4 (Rev. 05-01)

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002760
06/23/2011

09:59:07

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:59:50
Diagnostic Test	Pass	10:00:18
Air Blank	0.000	10:00:45
0.100 Cal Check	0.102	10:01:04
Air Blank	0.000	10:01:34
Subject Test	0.000	10:01:48
Air Blank	0.000	10:02:15
Five Minute Wait	ABT*	
Air Blank	0.000	10:03:47

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002760
06/23/2011

10:03:52

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	10:04:37
Air Blank	IPS**	10:04:45

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002760
06/23/2011

10:05:28

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:06:18
Diagnostic Test	Pass	10:06:45
Air Blank	0.000	10:07:13
0.100 Cal Check	0.101	10:07:32
Air Blank	0.000	10:08:01
Subject Test	0.031*	10:11:10
Air Blank	0.000	10:11:42
Air Blank	0.000	10:12:09
Subject Test	XXX**	10:12:32
Air Blank	0.000	10:13:05

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002760
06/23/2011

10:13:38

5 of 9 sc

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

BH

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:14:05
Cal Check	0.101	10:14:24
Air Blank	0.000	10:14:54

Sarah Cione
Operator's Signature

EXHIBIT G-2

6 of 9_{sc}

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TS4
6.27.2011

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 6/23/11 TIME 957
Intoxilyzer Serial # 80-002760 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: ARIZONA

Serial Number: 80-002760

06/23/2011

09:57:33

Standard Lot#: 20910100A1

Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:58:01
Cal Check	0.101	09:58:20
Air Blank	0.000	09:58:49

Sarah Cione

Operator's Signature

7 of 9 SC

154



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

8 of 9_{sc}

INTOXILYZER SERVICE RECORD
FOR
CALIBRATION GAS TANK CHANGE

154
6-27-2011

Date: 6/23/11

QAS: Sarah Cione

Intoxilyzer Serial Number: 80-002760

New Lot Number: 20910100A1

New Expiration Date: 8/1/2012

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."
Select "D" for dry gas and enter the following information:
Enter Target Value: 0.100
Enter Control Lot #: From the tank label
Enter Exp Date: From the tank label (mm/dd/yyyy)
Edit info Y/N? Y or N
Update Std Info? Y ENTER
"Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "Control Testing."
Select "D" for dry control test. After selection is made press
ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 100 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Sarah Cione


INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002760
06/23/2011

10:18:13

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:18:41
Cal Check	0.100	10:18:59
Air Blank	0.000	10:19:29

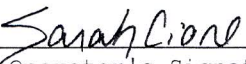

Operator's Signature