



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

B4
8-22-2011

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002760 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 8-18-11 Cal 1 0.100 AC Cal 2 0.100 AC

Between the dates of 7-21-11 and 8-18-11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.


Clare Conley
Clare L. Conley

8-18-11

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 18th, day of August, 2011.

Susan A. Hitchye
Notary Public

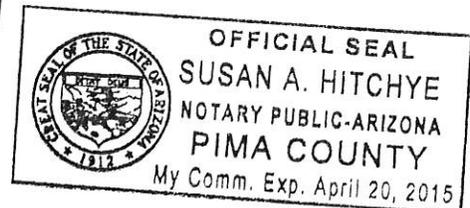

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1567
cm
TS4
8-22-2011

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 8-18-11 TIME 0851

INTOXILYZER SERIAL # SD-002760 LOCATION ODW DMU HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.

(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.

(✓) 2. Error recognition logic system functioning.

Not a Successfully Completed Test Sequence printed

(✓) 3. Proper sample recognition system.

Not a Successfully Completed Test Sequence printed

Deficient sample printed.

(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES NO

COMMENTS None.

SIGNATURE Clare Conley

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002760
08/18/2011

08:55:22

287
cu
B4

Standard Lot#: 04911100A1
Last Changed By: C. MORAND #24219

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:56:13
Diagnostic Test	Pass	08:56:41
Air Blank	0.000	08:57:08
0.100 Cal Check	0.101	08:57:27
Air Blank	0.000	08:57:57
Subject Test	0.000	08:58:23
Air Blank	ABT*	08:58:30
Air Blank	0.000	08:58:58

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002760
08/18/2011

08:59:18

357
CW

Standard Lot#: 04911100A1
Last Changed By: C. MORAND #24219

BY

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:00:03
Diagnostic Test	Pass	09:00:30
Air Blank	RFI*	09:00:37
Air Blank	IPS**	09:00:46

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002760
08/18/2011 09:01:13

487
CU
TS4

Standard Lot#: 04911100A1
Last Changed By: C. MORAND #24219

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST,SQA
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:01:59
Diagnostic Test	Pass	09:02:26
Air Blank	0.000	09:02:53
0.100 Cal Check	0.101	09:03:12
Air Blank	0.000	09:03:42
Subject Test	0.078*	09:06:50
Air Blank	0.000	09:07:23
Air Blank	0.000	09:07:50
Subject Test	XXX**	09:08:13
Air Blank	0.000	09:08:47

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002760
08/18/2011 09:09:06

507
CH

Standard Lot#: 04911100A1
Last Changed By: C. MORAND #24219

TS4

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:09:34
Cal Check	0.100	09:09:52
Air Blank	0.000	09:10:22

C. Conley

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

b.67
lu

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

B4
8.22.2011

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 8.18.11 TIME 0851
Intoxilyzer Serial # 80-002760 Location ODW DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002760
08/18/2011 08:51:26

7507
CW

Standard Lot#: 04911100A1
Last Changed By: C. MORAND #24219

TS4

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:51:54
Cal Check	0.100	08:52:13
Air Blank	0.000	08:52:42

Claude Conley

Operator's Signature



TUCSON POLICE DEPARTMENT
CRIME LABORATORY


72
Cu

INTOXILYZER SERVICE RECORD
FOR
CALIBRATION GAS TANK CHANGE

DM 8.22.2011

Date: 07-22-11
 QAS: V. Hernandez #42882 + C. Morand #24219
 Intoxilyzer Serial Number: 80-002760
 New Lot Number: 04911100A1
 New Expiration Date: 03-01-2013

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."
 Select "D" for dry gas and enter the following information:
 Enter Target Value: 0.100
 Enter Control Lot #: From the tank label
 Enter Exp Date: From the tank label (mm/dd/yyyy)
 Edit info Y/N? Y or N
 Update Std Info? Y ENTER
 "Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."
 Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 098 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove and sign printed record. Attach the record to the completed checklist.

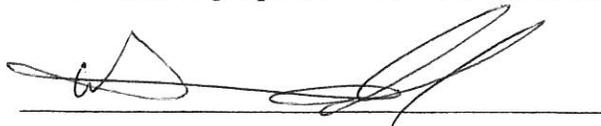SIGNATURE 

282
cu


154

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002760
07/22/2011 19:41:25

Standard Lot#: 04911100A1
Last Changed By: C. MORAND #24219

QAS: C. MORAND #24219
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	19:41:52
Cal Check	0.098	19:42:11
Air Blank	0.000	19:42:40

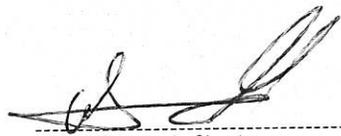

Operator's Signature