



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TS4
11-7-2011

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002760 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 10/14/11 Cal 1 0.099 AC Cal 2 0.098 AC

Between the dates of 9/16/11 and 10/14/11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -
The instrument was returned to service on -

Sarah Cione

Sarah L. Cione


10/27/11

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 27, day of October, 20 11.

A.M. Hernandez

Notary Public

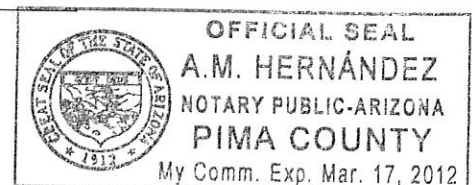

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

TS4
11.7.2011

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 10/14/11 TIME 1040

INTOXILYZER SERIAL # 80-002760 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.

(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0.000 AC.

(✓) 2. Error recognition logic system functioning.

Not a Successfully Completed Test Sequence printed

(✓) 3. Proper sample recognition system.

Not a Successfully Completed Test Sequence printed

Deficient sample printed.

(✓) 4. Standard Calibration Check standard 0.100 AC

Result: 0.098 AC

Instrument is operating properly and accurately. YES ✓ NO

COMMENTS

SIGNATURE

Sarah Cione

2 of 7sc

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002760
10/14/2011

10:45:12

TS4

Standard Lot#: 04911100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:45:57
Diagnostic Test	Pass	10:46:25
Air Blank	0.000	10:46:52
0.100 Cal Check	0.099	10:47:11
Air Blank	0.000	10:47:41
Subject Test	0.000	10:48:10
Air Blank	0.000	10:48:37
Five Minute Wait	ABT*	
Air Blank	0.000	10:49:12

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 7sc

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002760
10/14/2011

10:50:47

TS4

Standard Lot#: 04911100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	10:51:34
Air Blank	IPS**	10:51:44

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

4 of 7sc

154

INTOXILYZER 8000
 Location: ARIZONA
 Serial Number : 80-002760
 10/14/2011 10:52:23

Standard Lot#: 04911100A1
 Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
 TUCSON PD

Operator: S. CIONE #50616
 TUCSON PD

Subject: SQAP, SQAP
 DOB: 11/11/1111
 Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:53:06
Diagnostic Test	Pass	10:53:34
Air Blank	0.000	10:54:01
0.100 Cal Check	0.099	10:54:20
Air Blank	0.000	10:54:50
Subject Test	0.167*	10:57:58
Air Blank	0.000	10:58:31
Air Blank	0.000	10:58:58
Subject Test	XXX**	10:59:19
Air Blank	0.000	10:59:53

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: ARIZONA

Serial Number: 80-002760

10/14/2011

11:00:11

Standard Lot#: 04911100A1

Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:00:39
Cal Check	0.098	11:00:57
Air Blank	0.000	11:01:27

Sarah Cione

Operator's Signature

5 of 7 sc

TS4

EXHIBIT G-2

6 of 7 sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

BH
11.7.2011

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 10/14/11 TIME 1040
Intoxilyzer Serial # 80-002760 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002760
10/14/2011 10:40:53

Standard Lot#: 04911100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:41:21
Cal Check	0.099	10:41:39
Air Blank	0.000	10:42:09

Sarah Cione

Operator's Signature