



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 1306 W. Miracle Mile, Tucson, AZ 85705

TSY
 12-12-2011

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002760 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 12-8-11 Cal 1 0.097 AC Cal 2 0.097 AC

Between the dates of 11-14-11 and 12-8-11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
 The instrument was returned to service on N/A.

Clare Conley

Clare L. Conley


12-8-11

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 8th day of December, 2011.

Ann L. Beecroft

Notary Public

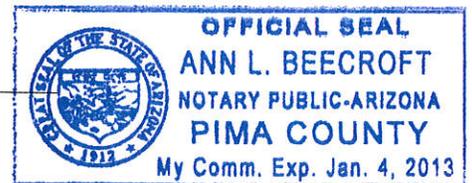

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

107
lu

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TS4
12-12-2011

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 12-8-11 TIME 0743

INTOXILYZER SERIAL # 80-002760 LOCATION ODW DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 097 AC

Instrument is operating properly and accurately. YES NO

COMMENTS None.

SIGNATURE Clare Conley

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002760
12/08/2011

08:13:03

2A7
cu

Standard Lot#: 04911100A1
Last Changed By: S. CIONE #50616

B4

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:13:53
Diagnostic Test	Pass	08:14:20
Air Blank	0.000	08:14:47
0.100 Cal Check	0.098	08:15:06
Air Blank	0.000	08:15:36
Subject Test	0.000	08:17:22
Air Blank	ABT*	08:17:28
Air Blank	IPS**	08:17:36

*Sequence Aborted
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002760
12/08/2011

08:17:41

Standard Lot#: 04911100A1
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	08:18:11
Air Blank	0.000	08:18:39

*RFI Detect

Not a Successfully Completed Test Sequence

3567
CW

PH

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002760
12/08/2011

08:19:01

Standard Lot#: 04911100A1
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:19:44
Diagnostic Test	Pass	08:20:11
Air Blank	0.000	08:20:38
0.100 Cal Check	0.098	08:20:57
Air Blank	0.000	08:21:27
Subject Test	0.074*	08:24:35
Air Blank	0.000	08:25:07
Air Blank	0.000	08:25:35
Subject Test	XXX**	08:26:59
Air Blank	0.000	08:27:33

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002760
12/08/2011 08:27:52

5007
en

Standard Lot#: 04911100A1
Last Changed By: S. CIONE #50616

T24

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:28:19
Cal Check	0.097	08:28:38
Air Blank	0.000	08:29:08

C. Conley

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

657
Cu

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

134
12.12.2011

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 12-8-11 TIME 0743
Intoxilyzer Serial # 80-002760 Location ODW DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 097 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002760
12/08/2011 07:43:00

787
lu

Standard Lot#: 04911100A1
Last Changed By: S. CIONE #50616

TS4

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	07:43:27
Cal Check	0.097	07:43:46
Air Blank	0.000	07:44:15

C. Conley

Operator's Signature