



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

TS4  
1-18-2012

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002760 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 1/5/12 Cal 1 0.098 AC Cal 2 0.097 AC

Between the dates of 12/8/11 and 1/5/12 the instrument was operating properly and accurately.


The instrument was taken out of service for repair on -.  
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione

1/11/12


ADPS QAS Permit # 33355

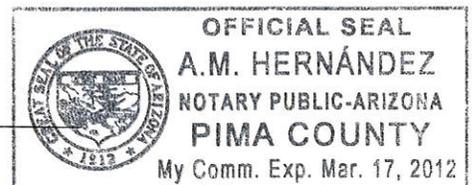

SUBSCRIBED AND SWORN to before me this 11, day of January, 20 12.

A.M. Hernandez

Notary Public

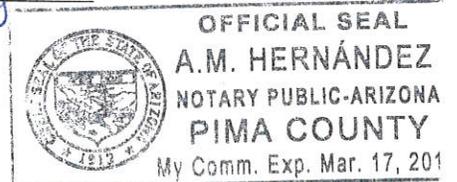

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

TSG  
1-18-2012

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 1/5/12 TIME 1418

INTOXILYZER SERIAL # 80-002760 LOCATION DUI HQ

(  ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (  ) 1. Clock time check.
- (  ) 2. Date check.

OPERATIONAL TESTS

- (  ) 1. Alcohol-free subject test result 0. 000 AC.
- (  ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (  ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (  ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 097 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department

2 of 7sc

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002760  
01/05/2012

14:20:25

TS4

Standard Lot#: 04911100A1  
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:21:08
Diagnostic Test	Pass	14:21:35
Air Blank	0.000	14:22:03
0.100 Cal Check	0.099	14:22:21
Air Blank	0.000	14:22:51
Subject Test	0.000	14:23:11
Air Blank	0.000	14:23:38
Five Minute Wait	ABT*	
Air Blank	0.000	14:24:46

\*Sequence Aborted

### Not a Successfully Completed Test Sequence

3 of 7<sub>sc</sub>

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002760  
01/05/2012 14:25:33

TS4

---

Standard Lot#: 04911100A1  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	RFI*	14:26:01
Air Blank	IPS**	14:26:08

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

4 of 7 SC

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002760  
01/05/2012

14:27:27

TSh

Standard Lot#: 04911100A1  
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:28:21
Diagnostic Test	Pass	14:28:48
Air Blank	0.000	14:29:16
0.100 Cal Check	0.098	14:29:34
Air Blank	0.000	14:30:04
Subject Test	0.133*	14:33:12
Air Blank	0.000	14:33:46
Air Blank	0.000	14:34:13
Subject Test	XXX**	14:34:31
Air Blank	0.000	14:35:05

\*Deficient Sample  
\*\*Invalid Sample  
Begin new deprivation period

**Not a Successfully Completed Test Sequence**

5 of 7sc

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002760  
01/05/2012 14:35:42

ASH

-----  
Standard Lot#: 04911100A1  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	14:36:09
Cal Check	0.097	14:36:28
Air Blank	0.000	14:36:58

*Sarah Cione*

\_\_\_\_\_  
Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD CALIBRATION CHECK PROCEDURE

TS4  
1-18-2012

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 1/5/12 TIME 1418  
Intoxilyzer Serial # 80-002760 Location DUI HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 098 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

7 of 7sc

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002760  
01/05/2012 14:18:13

TS4

-----  
Standard Lot#: 04911100A1  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	14:18:40
Cal Check	0.098	14:18:59
Air Blank	0.000	14:19:28

Sarah Cione  
Operator's Signature