



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 1306 W. Miracle Mile, Tucson, AZ 85705

TS4
 1-28-2013

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002760 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 1/3/13 Cal 1 0.100 AC Cal 2 0.100 AC

Between the dates of 12/6/12 and 1/3/13 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
 The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


1/7/13

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 7th, day of January, 2013.

Elma Escarcega

Notary Public

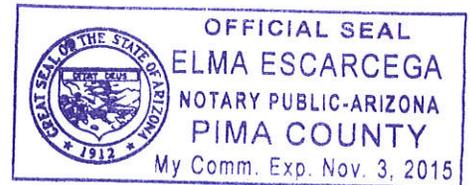

*Includes documents specified in A.A.C. R13-10-104.A.5

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 8_{sc}
By
1-28-2013

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 1/3/13 TIME 1136

INTOXILYZER SERIAL # 80-002760 LOCATION DUI HQ

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

() 1. Clock time check.

() 2. Date check.

OPERATIONAL TESTS

() 1. Alcohol-free subject test result 0. 000 AC.

() 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed

() 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.

() 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002760
Core Version: 8105.48
01/03/2013

11:38:36

2 of 8_{sc}
Bh

Standard Lot#: 23511100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1980

Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:39:27
Diagnostic Test	Pass	11:40:01
Air Blank	0.000	11:40:29
0.100 Cal Check	0.100	11:40:50
Air Blank	0.000	11:41:20
Subject Test	0.000	11:41:53
Air Blank	0.000	11:42:21
Wait	ABT*	
Air Blank	0.000	11:43:19

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002760
Core Version: 8105.48
01/03/2013

11:43:46

3 of 8 sc
TBH

Standard Lot#: 23511100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:44:35
Diagnostic Test	Pass	11:45:10
Air Blank	0.000	11:45:37
0.100 Cal Check	0.100	11:45:58
Air Blank	0.000	11:46:28
Subject Test	RFI*	11:46:42
Air Blank	IPS**	11:46:50

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

4 of 8 sc
TS

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002760
Core Version: 8105.48
01/03/2013

11:59:29

Standard Lot#: 23511100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:00:12
Diagnostic Test	Pass	12:00:47
Air Blank	0.000	12:01:15
0.100 Cal Check	0.101	12:01:36
Air Blank	0.000	12:02:05
Subject Test	XXX*	12:02:34
Air Blank	0.000	12:03:10

*Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

5 of 8 sc
TS4

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002760
Core Version: 8105.48
01/03/2013

12:03:33

Standard Lot#: 23511100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:04:20
Diagnostic Test	Pass	12:04:55
Air Blank	0.000	12:05:22
0.100 Cal Check	0.101	12:05:43
Air Blank	0.000	12:06:12
Subject Test	0.201*	12:09:20
Air Blank	0.000	12:09:53
Air Blank	0.000	12:10:21
Subject Test	ABT**	12:10:36
Air Blank	0.000	12:11:04

*Deficient Sample
**Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002760
Core Version: 8105.48
01/03/2013

12:11:44

6 of 85
TS4

Standard Lot#: 23511100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	12:12:12
0.100 Cal Check	0.100	12:12:33
Air Blank	0.000	12:13:03

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

7 of 8 sc
TSK
1.28.2013

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 1/3/13 TIME 1136
Intoxilyzer Serial # 80-002760 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002760
Core Version: 8105.48
01/03/2013

11:36:29

8 of 8 sc
TS4

Standard Lot#: 23511100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:36:58
0.100 Cal Check	0.100	11:37:19
Air Blank	0.000	11:37:48

Preventative Maintenance
Performed -- Timer Reset