



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 1306 W. Miracle Mile, Tucson, AZ 85705

TS4
 2-19-2013

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002760 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 1/31/13 Cal 1 0.099 AC Cal 2 0.099 AC

Between the dates of 1/3/13 and 1/31/13 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
 The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


2/14/13

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 14th, day of February, 2013.

Elma Escarcega

Notary Public

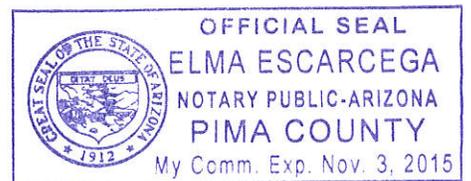

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 7_{sc}
TB4
2-19-2013

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 1/31/13 TIME 1049

INTOXILYZER SERIAL # 80-002760 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.

(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.

(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed

(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.

(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

2 of 7_{sc}
TS4

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002760
Core Version: 8105.48
01/31/2013 10:56:09

Standard Lot#: 23511100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:57:06
Diagnostic Test	Pass	10:57:41
Air Blank	0.000	10:58:08
0.100 Cal Check	0.100	10:58:29
Air Blank	0.000	10:58:59
Subject Test	0.000	10:59:18
Air Blank	0.000	10:59:45
Wait	ABT*	
Air Blank	0.000	11:01:08

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 7_{sc}
B4

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002760
Core Version: 8105.48
01/31/2013 11:05:05

Standard Lot#: 23511100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:05:49
Diagnostic Test	Pass	11:06:24
Air Blank	0.000	11:06:51
0.100 Cal Check	0.100	11:07:12
Air Blank	0.000	11:07:41
Subject Test	RFI*	11:07:53
Air Blank	IPS**	11:08:08

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

4 of 7_{sc}
B4

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002760
Core Version: 8105.48
01/31/2013

11:11:52

Standard Lot#: 23511100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:12:49
Diagnostic Test	Pass	11:13:24
Air Blank	0.000	11:13:52
0.100 Cal Check	0.101	11:14:13
Air Blank	0.000	11:14:43
Subject Test	0.151*	11:17:50
Air Blank	0.000	11:18:23
Air Blank	0.000	11:18:50
Subject Test	XXX**	11:19:15
Air Blank	0.000	11:19:51

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002760
Core Version: 8105.48
01/31/2013 11:20:29

5 of 7_{sc}
B4

Standard Lot#: 23511100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:20:58
0.100 Cal Check	0.099	11:21:19
Air Blank	0.000	11:21:49

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6 of 7 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TSh
2.19.2013

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 1/31/13 TIME 1049
Intoxilyzer Serial # 80-002760 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002760
Core Version: 8105.48
01/31/2013 10:49:29

7 of 7sc

TS4

Standard Lot#: 23511100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:49:57
0.100 Cal Check	0.099	10:50:18
Air Blank	0.000	10:50:48

Preventative Maintenance
Performed -- Timer Reset