



CERTIFICATION by CUSTODIAN of RECORDS  
 Tucson Police Department Crime Laboratory  
 1306 W. Miracle Mile, Tucson, AZ 85705

TSy  
 1.14.2014

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002760 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 1/2/14 Cal 1 0.097 AC Cal 2 0.096 AC

Between the dates of 12/5/13 and 1/2/14 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on 1/2/14.  
 The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


1/9/14

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 9<sup>th</sup>, day of January, 2014.

Elma Escarcega

Notary Public

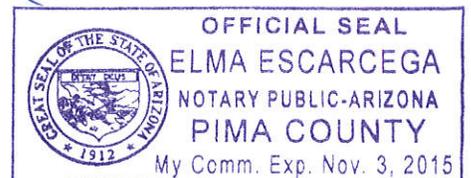

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 9  
SC  
TS4  
1-14-2014

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 1/2/14 TIME 926

INTOXILYZER SERIAL # 80-002760 LOCATION DUI HQ

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 096 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE Sarah Cione

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002760  
Core Version: 8105.48  
01/02/2014

09:30:07

2 of 9  
SC  
B4

---

Standard Lot#: 30412100A1  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	09:30:50
Diagnostic Test	Pass	09:31:25
Air Blank	0.000	09:31:52
0.100 Cal Check	0.096	09:32:13
Air Blank	0.000	09:32:43
Subject Test	0.000	09:33:33
Air Blank	0.000	09:34:01
Wait	ABT*	
Air Blank	0.000	09:34:51

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002760  
Core Version: 8105.48  
01/02/2014

09:35:33

3 of 9  
SC

TSh

---

Standard Lot#: 30412100A1  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	09:36:35
Diagnostic Test	Pass	09:37:10
Air Blank	0.000	09:37:37
0.100 Cal Check	0.097	09:37:58
Air Blank	0.000	09:38:28
Subject Test	RFI*	09:39:01
Air Blank	IPS**	09:39:12

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002760  
Core Version: 8105.48  
01/02/2014

09:41:55

4 of 9  
SC  
TS4

---

Standard Lot#: 30412100A1  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	09:42:41
Diagnostic Test	Pass	09:43:16
Air Blank	0.000	09:43:43
0.100 Cal Check	0.097	09:44:04
Air Blank	0.000	09:44:34
Subject Test	0.365*	09:47:41
Air Blank	0.000	09:48:14
Air Blank	0.000	09:48:42
Subject Test	0.612	09:49:01
Air Blank	0.000	09:49:42
Wait		
Air Blank	0.000	09:54:29
Subject Test	XXX**	09:55:45
Air Blank	0.000	09:56:18

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002760  
Core Version: 8105.48  
01/02/2014 09:57:19

5 of 9  
SC  
TSh

-----  
Standard Lot#: 30412100A1  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:57:47
0.100 Cal Check	0.096	09:58:08
Air Blank	0.000	09:58:38

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 9  
SC  
TS4

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 1/2/14 TIME 926  
Intoxilyzer Serial # 80-002760 Location DUI HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0.100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0.097 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002760  
Core Version: 8105.48  
01/02/2014 09:26:40

7 of 9  
SC  
TSh

-----  
Standard Lot#: 30412100A1  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:27:08
0.100 Cal Check	0.097	09:27:29
Air Blank	0.000	09:27:59

Preventative Maintenance  
Performed -- Timer Reset



8 of 9 sc  
B4

### Service Evaluation Form

This form MUST be completed and enclosed with instrument to be serviced.  
Failure to complete and return this form may cause delays in service.

(Note: please ship items in their original shipping container(s) or a similar protective box.)

1. Contact information: Customer Number 105075 (contact Customer Service)  
 Name Terry Gallegos Phone: (520) 837-7677  
 Fax: (520) 791-4642 Email: Terry.Gallegos@tucsonaz.gov

2. Bill to Address: Ship to Address:  
Finance / Accounts Payable Tucson Police Department Crime Lab  
Po Box 27450 1306 W. Miracle Mile  
Tucson, AZ 85726-7450 Tucson, AZ 85705

3. Serial Number: 80-002760 Instrument Model: I8000

4. Detailed Description of Problem:  
Measured past 0.600 - did not give "range exceeded" message

\*\*\*Hazardous Material Warning! - DO NOT return gas cylinders with instrument!\*\*\*

I Authorize Repairs Up To:  All  \$250  \$500  \$750  Other \$ \_\_\_\_\_  
 Purchase Order Number (attach a copy of P.O. if applicable) 34770

Authorized By:

Terry Gallegos Toxicology Coordinator  
 Name (Please Print) Title  
[Signature] 1-14-2014  
 Signature Date

Ship to:  
 CMI, Inc.  
 Attn: Service Dept.  
 316 East Ninth Street  
 Owensboro, KY 42303

No, please send estimate before repairs are made.  
 Note: An estimate will be faxed before performing any repairs and may cause delays in service.  
 An evaluation fee (\$79.00 or actual costs) will apply to estimates that are not repaired.



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

9 of 9 sc  
TS4

INTOXILYZER SERVICE RECORD

Date: 1/2/14

QAS: Sarah Cione

Intoxilyzer Serial Number: 80-002760

Reason For Service: Registered past 0.600 - did not give range  
exceeded message

Date Out of Service: 1/2/14

Calibration Check completed prior to removal from service

Date: 1/2/14

SQAP completed prior to removal from service

Date: 1/2/14

Corrective Action:

Sent to CMI, Inc.

Data downloaded prior to sending in to CMI