



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TS4
8.4.2014

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002760 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 7/18/14 Cal 1 0.097 AC Cal 2 0.097 AC

Between the dates of 6/18/14 and 7/18/14 the instrument was operating properly and accurately.


The instrument was taken out of service for repair on -.
The instrument was returned to service on -.

Sarah Cione
Sarah L. Cione

7/21/14
ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 21st, day of July, 2014.

Elma Escarcega
Notary Public

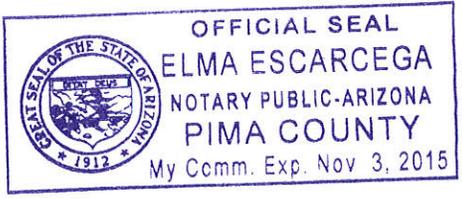

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7
JC
B4
8.4.2014

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 7/18/14 TIME 959

INTOXILYZER SERIAL # 80-002760 LOCATION Crime Lab

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 097 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002760
Core Version: 8105.48
07/18/2014

10:01:37

2 of 7
sc
B4

Standard Lot#: 19613100A4
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 10:02:40 |
| Diagnostic Test | Pass | 10:03:15 |
| Air Blank | 0.000 | 10:03:43 |
| 0.100 Cal Check | 0.097 | 10:04:04 |
| Air Blank | 0.000 | 10:04:33 |
| Subject Test | 0.000 | 10:04:55 |
| Air Blank | 0.000 | 10:05:22 |
| Wait | ABT* | |
| Air Blank | 0.000 | 10:05:56 |

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002760
Core Version: 8105.48
07/18/2014

10:06:03

3 of 7 SC
BH

Standard Lot#: 19613100A4
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 10:07:05 |
| Diagnostic Test | Pass | 10:07:40 |
| Air Blank | 0.000 | 10:08:07 |
| 0.100 Cal Check | 0.098 | 10:08:28 |
| Air Blank | 0.000 | 10:08:58 |
| Subject Test | RFI* | 10:09:21 |
| Air Blank | IPS** | 10:09:29 |

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002760
Core Version: 8105.48
07/18/2014

10:09:36

4 of 7
sc
B4

Standard Lot#: 19613100A4
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 10:10:36 |
| Diagnostic Test | Pass | 10:11:11 |
| Air Blank | 0.000 | 10:11:39 |
| 0.100 Cal Check | 0.097 | 10:11:59 |
| Air Blank | 0.000 | 10:12:29 |
| Subject Test | 0.177* | 10:15:36 |
| Air Blank | 0.000 | 10:16:08 |
| Air Blank | 0.000 | 10:16:36 |
| Subject Test | XXX** | 10:17:26 |
| Air Blank | 0.000 | 10:17:59 |

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002760
Core Version: 8105.48
07/18/2014

10:18:37

5 of 7
sc
TS4

Standard Lot#: 19613100A4
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 10:19:06 |
| 0.100 Cal Check | 0.097 | 10:19:27 |
| Air Blank | 0.000 | 10:19:56 |

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7 sc
TS4
8.4.2014

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 7/18/14 TIME 959
Intoxilyzer Serial # 80-002760 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 097 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002760
Core Version: 8105.48
07/18/2014 09:59:01

7 of 7 sc
B4
8.4.2014

Standard Lot#: 19613100A4
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 09:59:30 |
| 0.100 Cal Check | 0.097 | 09:59:51 |
| Air Blank | 0.000 | 10:00:20 |

Preventative Maintenance
Performed -- Timer Reset