



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

TSh  
9-5-2014

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002760 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 8/15/14 Cal 1 0.096 AC Cal 2 0.097 AC

Between the dates of 7/18/14 and 8/15/14 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


8/19/14

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 19th, day of August, 20 14.

Elma Escarcega

Notary Public

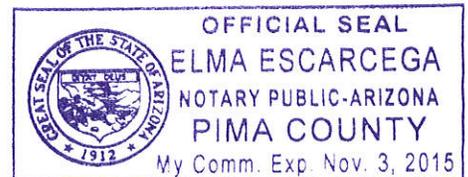

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1 of 7  
sc  
BY  
9.5.2014

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 8/15/14 TIME 1517

INTOXILYZER SERIAL # 80-002760 LOCATION Crime Lab

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 097 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Sarah Cione

2 of 7  
sc  
bh

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002760  
Core Version: 8105.48  
08/15/2014 15:19:05

-----  
Standard Lot#: 19613100A4  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Subject: SOAP, SOAP  
DOB: 11/11/1980  
Sex: F Weight: 111  
-----

15 Minute Deprivation Period? Yes

-----

Test	g/210L	Time
Air Blank	0.000	15:19:47
Diagnostic Test	Pass	15:20:24
Air Blank	0.000	15:20:51
0.100 Cal Check	0.097	15:21:12
Air Blank	0.000	15:21:41
Subject Test	0.000	15:22:02
Air Blank	0.000	15:22:30
Wait	ABT*	
Air Blank	0.000	15:23:10

\*Sequence Aborted

Not a Successfully  
Completed  
Test Sequence

3 of 7  
sc  
Rsh

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002760  
Core Version: 8105.48  
08/15/2014 15:24:10

-----  
Standard Lot#: 19613100A4  
Last Changed By: S. CIONE #50616  
-----

DAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111  
-----

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:24:56
Diagnostic Test	Pass	15:25:33
Air Blank	0.000	15:26:01
0.100 Cal Check	0.097	15:26:21
Air Blank	0.000	15:26:51
Subject Test	RFI*	15:27:59
Air Blank	IPS**	15:28:07

\*RFI Detect  
\*\*Improper Sample

Not a Successfully  
Completed  
Test Sequence

4 of 7 sc  
Bh

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002760  
Core Version: 8105.48  
08/15/2014 15:29:05

-----  
Standard Lot#: 19613100A4  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Subject: SOAP, SOAP  
DOB: 11/11/1980  
Sex: F Weight: 111  
-----

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:29:52
Diagnostic Test	Pass	15:30:30
Air Blank	0.000	15:30:57
0.100 Cal Check	0.097	15:31:18
Air Blank	0.000	15:31:47
Subject Test	0.313*	15:34:55
Air Blank	0.000	15:35:27
Air Blank	0.000	15:35:54
Subject Test	XXX**	15:36:10
Air Blank	0.000	15:36:43

\*Deficient Sample  
\*\*Invalid Sample  
Begin new deprivation period

Not a Successfully  
Completed  
Test Sequence

5 of 7  
SC  
Bh

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002760  
Core Version: 8105.48  
08/15/2014 15:37:48

-----  
Standard Lot#: 19613100A4  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	15:38:16
0.100 Cal Check	0.097	15:38:37
Air Blank	0.000	15:39:07

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7  
sc  
By  
9.5.2014

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 8/15/14 TIME 1517  
Intoxilyzer Serial # 80-002760 Location Crime Lab

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 096 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

7 of 7 sc  
TS4

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002760  
Core Version: 8105.48  
08/15/2014 15:17:09

-----  
Standard Lot#: 19613100A4  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	15:17:37
0.100 Cal Check	0.096	15:17:58
Air Blank	0.000	15:18:27

Preventative Maintenance  
Performed -- Timer Reset