



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

SC  
12/19/14

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002760 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 12-4-2014 Cal 1 0.097 AC Cal 2 0.096 AC

Between the dates of 11-6-2014 and 12-4-2014 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.  
The instrument was returned to service on —.


Terry S. Gallegos  
Terry S. Gallegos

12-17-2014

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 17<sup>th</sup>, day of December, 2014.

Elma Escarcega  
Notary Public

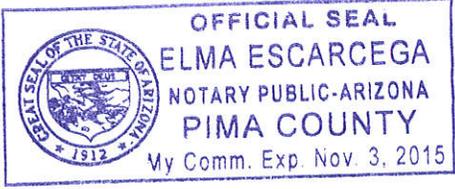

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

10 of 7  
TS4

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

SC  
12/19/14

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 12-4-2014 TIME 1554  
INTOXILYZER SERIAL # 80-002760 LOCATION ODW-DUI-HQ

(  ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (  ) 1. Clock time check.
- (  ) 2. Date check.

OPERATIONAL TESTS

- (  ) 1. Alcohol-free subject test result 0. 000 AC.
- (  ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (  ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (  ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 096 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE

Terry Gallegos

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002760  
Core Version: 8105.48  
12/04/2014

15:57:44

2 of 7  
TS6  
5C

Standard Lot#: 19613100A4  
Last Changed By: S. CIONE #50616

QAS: T. GALLEGOS #28847  
TUCSON PD

Operator: T. GALLEGOS #28847  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:58:58
Diagnostic Test	Pass	15:59:33
Air Blank	0.000	16:00:00
0.100 Cal Check	0.097	16:00:21
Air Blank	0.000	16:00:51
Subject Test	0.000	16:01:22
Air Blank	0.000	16:01:49
Wait	ABT*	
Air Blank	0.000	16:02:23

\*Sequence Aborted

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002760  
Core Version: 8105.48  
12/04/2014

16:02:52

30F7  
T54

5C

Standard Lot#: 19613100A4  
Last Changed By: S. CIONE #50616

QAS: T. GALLEGOS #28847  
TUCSON PD

Operator: T. GALLEGOS #28847  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:04:14
Diagnostic Test	Pass	16:04:49
Air Blank	0.000	16:05:17
0.100 Cal Check	RFI*	16:05:33
Air Blank	IPS**	16:05:45

\*RFI Detect  
\*\*Improper Sample

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002760  
Core Version: 8105.48  
12/04/2014

16:06:27

40F7  
TS4

5C

Standard Lot#: 19613100A4  
Last Changed By: S. CIONE #50616

QAS: T. GALLEGOS #28847  
TUCSON PD

Operator: T. GALLEGOS #28847  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:07:17
Diagnostic Test	Pass	16:07:52
Air Blank	0.000	16:08:19
0.100 Cal Check	0.097	16:08:40
Air Blank	0.000	16:09:10
Subject Test	0.072*	16:12:17
Air Blank	0.000	16:12:50
Air Blank	0.000	16:13:17
Subject Test	XXX**	16:13:32
Air Blank	0.000	16:14:06

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002760  
Core Version: 8105.48  
12/04/2014 16:14:37

507  
754

5c

-----  
Standard Lot#: 19613100A4  
Last Changed By: S. CIONE #50616  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	$\sigma/210L$	Time
Air Blank	0.000	16:15:06
0.100 Cal Check	0.096	16:15:27
Air Blank	0.000	16:15:57

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

6 of 7  
TS4

ARIZONA DEPARTMENT OF PUBLIC SAFETY

SC

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

12/19/14

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 12-4-2014 TIME 1554  
Intoxilyzer Serial # 80-002760 Location DDW-DUI-HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 097 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002760  
Core Version: 8105.48  
12/04/2014 15:54:11

7 of 7  
Bh

-----  
Standard Lot#: 19613100A4  
Last Changed By: S. CIONE #50616  
-----

SC

-----  
QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

-----  
Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	15:54:41
0.100 Cal Check	0.097	15:55:02
Air Blank	0.000	15:55:32

Preventative Maintenance  
Performed -- Timer Reset