



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002760 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4  
7-6-2015

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 5/28/2015 Cal 1 0.099 AC Cal 2 0.097 AC

Between the dates of 4/30/2015 and 5/28/2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on                     .  
The instrument was returned to service on                     .

William L Loomis

William L. Loomis


6/8/2015

ADPS QAS Permit # 33720

SUBSCRIBED AND SWORN to before me this 8<sup>th</sup>, day of June, 2015.

Elma Escarcega

Notary Public

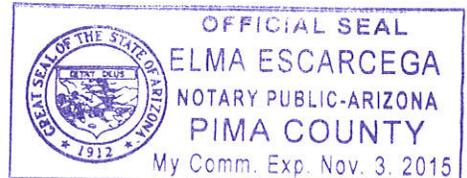

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7  
WU  
TS4 7-6-2015

QA SPECIALIST W. LOOMIS 53545 AGENCY Tucson Police Department

DATE 5/28/2015 TIME 14:53

INTOXILYZER SERIAL # 80-002760 LOCATION CRIME LAB

(  ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(  ) 1. Clock time check.

(  ) 2. Date check.

OPERATIONAL TESTS

(  ) 1. Alcohol-free subject test result 0. 000 AC.

(  ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed

(  ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.

(  ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 097 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Will Loomis

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002760  
Core Version: 8105.48  
05/28/2015 14:55:16

Standard Lot#: 25414100A4  
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545  
TUCSON PD

Operator: W. LOOMIS #53545  
TUCSON PD

Subject: SQAP,TEST  
DOB: 11/11/1987  
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:56:06
Diagnostic Test	Pass	14:56:41
Air Blank	0.000	14:57:08
0.100 Cal Check	0.099	14:57:29
Air Blank	0.000	14:57:59
Subject Test	0.000	14:58:16
Air Blank	0.000	14:58:43
Wait	ABT*	
Air Blank	0.000	14:59:18

\*Sequence Aborted

2 of 7  
WLL  
TSh

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002760  
Core Version: 8105.48  
05/28/2015 14:59:23

Standard Lot#: 25414100A4  
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545  
TUCSON PD

Operator: W. LOOMIS #53545  
TUCSON PD

Subject: SQAP,TEST  
DOB: 11/11/1987  
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:00:18
Diagnostic Test	Pass	15:00:53
Air Blank	0.000	15:01:20
0.100 Cal Check	RFI*	15:01:34
Air Blank	IPS**	15:01:42

\*RFI Detect  
\*\*Improper Sample

3 of 7  
WLL  
BH

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002760  
Core Version: 8105.48  
05/28/2015 15:01:48

4 of 7  
WLL  
TSG

---

Standard Lot#: 25414100A4  
Last Changed By: T. GALLEGOS #28847

---

QAS: W. LOOMIS #53545  
TUCSON PD

---

Operator: W. LOOMIS #53545  
TUCSON PD

---

Subject: SQAP,TEST  
DOB: 11/11/1987  
Sex: M Weight: 222

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	15:02:43
Diagnostic Test	Pass	15:03:18
Air Blank	0.000	15:03:45
0.100 Cal Check	0.099	15:04:06
Air Blank	0.000	15:04:36
Subject Test	0.146*	15:07:43
Air Blank	0.000	15:08:16
Air Blank	0.000	15:08:43
Subject Test	XXX**	15:09:05
Air Blank	0.000	15:09:39

\*Deficient Sample  
\*\*Invalid Sample  
Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002760  
Core Version: 8105.48  
05/28/2015 15:09:54

5 of 7  
WLL  
BH

-----  
Standard Lot#: 25414100A4  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: W. LOOMIS #53545  
TUCSON PD  
-----

Operator: W. LOOMIS #53545  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	15:10:23
0.100 Cal Check	0.097	15:10:44
Air Blank	0.000	15:11:13

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD CALIBRATION CHECK PROCEDURE

6 of 7  
will  
TS4  
7.6.2015

QA Specialist W. LOOMIS 53545 AGENCY Tucson Police Department  
DATE 5/28/2015 TIME 14:53  
Intoxilyzer Serial # 80-002760 Location CRIME LAB

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 099 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Will Loomis

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002760  
Core Version: 8105.48  
05/28/2015 14:53:14

7 of 7  
wll  
TS4

-----  
Standard Lot#: 25414100A4  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: W. LOOMIS #53545  
TUCSON PD  
-----

Operator: W. LOOMIS #53545  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	14:53:44
0.100 Cal Check	0.099	14:54:04
Air Blank	0.000	14:54:34

Preventative Maintenance  
Performed -- Timer Reset