



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 1306 W. Miracle Mile, Tucson, AZ 85705

TJG
 8-6-2015

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002760 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 7/24/15 Cal 1 0.097 AC Cal 2 0.095 AC

Between the dates of 6/25/15 and 7/24/15 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
 The instrument was returned to service on -.

Amy L. Autret
 Amy L. Autret


7/29/15

ADPS QAS Permit # 33718

SUBSCRIBED AND SWORN to before me this 29th, day of July, 2015.

Elma Escarcega

Notary Public

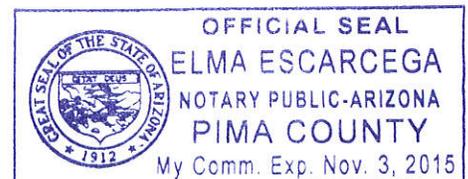

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1088
JLB
TS4
8.6.2015

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Amy Autret 53546 AGENCY Tucson Police Department

DATE 7/24/15 TIME 1029

INTOXILYZER SERIAL # 80-002760 LOCATION ODE - Intox room

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 095 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE

Amy Autret

2068
JLB
TS4

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002760
Core Version: 8105.48
07/24/2015

10:31:15

Standard Lot#: 25414100A4
Last Changed By: T. GALLEGOS #28847

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:32:12
Diagnostic Test	Pass	10:32:47
Air Blank	0.000	10:33:14
0.100 Cal Check	0.097	10:33:35
Air Blank	0.000	10:34:04
Subject Test	0.000	10:34:18
Air Blank	0.000	10:34:46
Wait	ABT*	
Air Blank	0.000	10:35:20

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 8
JL
134

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002760
Core Version: 8105.48
07/24/2015 10:35:28

Standard Lot#: 25414100A4
Last Changed By: T. GALLEGOS #28847

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:36:26
Diagnostic Test	Pass	10:37:01
Air Blank	0.000	10:37:28
0.100 Cal Check	RFI*	10:37:40
Air Blank	IPS**	10:37:48

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

4 of 8
JLB
TSh

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002760
Core Version: 8105.48
07/24/2015

10:37:59

Standard Lot#: 25414100A4
Last Changed By: T. GALLEGOS #28847

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:38:53
Diagnostic Test	Pass	10:39:28
Air Blank	0.000	10:39:56
0.100 Cal Check	0.098	10:40:17
Air Blank	0.000	10:40:47
Subject Test	0.314*	10:43:54
Air Blank	0.000	10:44:32
Air Blank	RFI**	10:44:48
Air Blank	0.000	10:45:16

*Deficient Sample
**RFI Detect

Several officers
came into room -
unexpected RFI -
from nearby fax machine
being used, cell phones,
or police radios?
Repeated Deficient Sample
& Invalid Sample Sequence.

Not a Successfully Completed Test Sequence

JLB

5 of 8
MB
TB4

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002760
Core Version: 8105.48
07/24/2015 10:46:21

Standard Lot#: 25414100A4
Last Changed By: T. GALLEGOS #28847

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:47:16
Diagnostic Test	Pass	10:47:51
Air Blank	0.000	10:48:19
0.100 Cal Check	0.096	10:48:39
Air Blank	0.000	10:49:10
Subject Test	0.276*	10:52:17
Air Blank	0.000	10:52:55
Air Blank	0.000	10:53:22
Subject Test	XXX**	10:53:51
Air Blank	0.000	10:54:29

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002760
Core Version: 8105.48
07/24/2015 10:54:53

leafs
Ad
Bh

Standard Lot#: 25414100A4
Last Changed By: T. GALLEGOS #28847

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:55:22
0.100 Cal Check	0.095	10:55:43
Air Blank	0.000	10:56:13

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

7068
[Signature]

TS4
8.6.2015

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Amy Antret 53546 AGENCY Tucson Police Department

DATE 7/24/15 TIME 1029

Intoxilyzer Serial # 80-002760 Location OPE-Intox room

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 097 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE *[Signature]*

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002760
Core Version: 8105.48
07/24/2015 10:29:45

8 of 8
dd

TBY

Standard Lot#: 25414100A4
Last Changed By: T. GALLEGOS #28847

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:30:14
0.100 Cal Check	0.097	10:30:35
Air Blank	0.000	10:31:04

Preventative Maintenance
Performed -- Timer Reset