



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

TS4  
11-3-2015

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002760 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 10/15/2015 Cal 1 0.097 AC Cal 2 0.095 AC

Between the dates of 9/17/2015 and 10/15/2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on \_\_\_\_\_.  
The instrument was returned to service on \_\_\_\_\_.

William L Loomis

William L. Loomis

10/22/2015


ADPS QAS Permit # 33720

SUBSCRIBED AND SWORN to before me this 22nd, day of October, 2015.

Elma Escarcega

Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

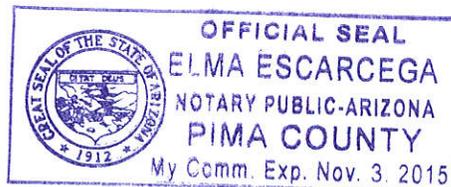

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7  
WLL  
TSK  
11-3-2015

QA SPECIALIST W. LOOMIS 53545 AGENCY Tucson Police Department  
DATE 10/15/2015 TIME 09:18  
INTOXILYZER SERIAL # 80-002760 LOCATION ODE - SGT'S OFFICE

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 095 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Will Loomis

2 of 7  
wll  
B4

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002760  
Core Version: 8105.48  
10/15/2015 09:21:17

-----  
Standard Lot#: 25414100A4  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: W. LOOMIS #53545  
TUCSON PD  
-----

Operator: W. LOOMIS #53545  
TUCSON PD  
-----

Subject: SOAP, TEST  
DOB: 11/11/1987  
Sex: M Weight: 222  
-----

15 Minute Deprivation Period? Yes  
-----

Test	g/210L	Time
Air Blank	0.000	09:22:09
Diagnostic Test	Pass	09:22:47
Air Blank	0.000	09:23:14
0.100 Cal Check	0.096	09:23:35
Air Blank	0.000	09:24:05
Subject Test	0.000	09:24:40
Air Blank	0.000	09:25:07
Wait	ABT*	
Air Blank	0.000	09:25:43

\*Sequence Aborted

Not a Successfully  
Completed  
Test Sequence

3 of 7  
wv  
B4

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002760  
Core Version: 8105.48  
10/15/2015 09:29:23

-----  
Standard Lot#: 25414100A4  
Last Changed By: T. GALLEGOS #28847  
-----

OAS: W. LOOMIS #53545  
TUCSON PD  
-----

Operator: W. LOOMIS #53545  
TUCSON PD  
-----

Subject: SQAP, TEST  
DOB: 11/11/1987  
Sex: M Weight: 222  
-----

15 Minute Deprivation Period? Yes  
-----

Test	g/210L	Time
Air Blank	0.000	09:30:19
Diagnostic Test	Pass	09:30:57
Air Blank	0.000	09:31:24
0.100 Cal Check	RFI*	09:31:38
Air Blank	IPS**	09:31:46

\*RFI Detect  
\*\*Improper Sample

Not a Successfully  
Completed  
Test Sequence

4 of 7  
WLL  
TSH

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002760  
Core Version: 8105.48  
10/15/2015 09:33:44

-----  
Standard Lot#: 25414100A4  
Last Changed By: T. GALLEGOS #28847  
-----

OAS: W. LOOMIS #53545  
TUCSON PD  
-----

Operator: W. LOOMIS #53545  
TUCSON PD  
-----

Subject: SOAP, TEST  
DOB: 11/11/1987  
Sex: M Weight: 222  
-----

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:34:39
Diagnostic Test	Pass	09:35:17
Air Blank	0.000	09:35:44
0.100 Cal Check	0.096	09:36:05
Air Blank	0.000	09:36:35
Subject Test	0.186*	09:39:42
Air Blank	0.000	09:40:15
Air Blank	0.000	09:40:42
Subject Test	XXX**	09:40:56
Air Blank	0.000	09:41:30

\*Deficient Sample  
\*\*Invalid Sample  
Begin new deprivation period

Not a Successfully  
Completed  
Test Sequence

5 of 7  
WW  
TBH

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002760  
Core Version: 8105.48  
10/15/2015 09:43:51

-----  
Standard Lot#: 25414100A4  
Last Changed By: T. GALLEGOS #28847  
-----

GAS: W. LOOMIS #53545  
TUCSON PD  
-----

Operator: W. LOOMIS #53545  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:44:21
0.100 Cal Check	0.095	09:44:42
Air Blank	0.000	09:45:11

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7  
WU  
TSY  
11.3.2015

QA Specialist W. LOOMIS 53545 AGENCY Tucson Police Department

DATE 10/15/2015 TIME 09:18

Intoxilyzer Serial # 80-002760 Location ODE - SGT'S OFFICE

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C.
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 097 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Will Loomis

7 of 7  
wll  
B4

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002760  
Core Version: 8105.48  
10/15/2015 09:18:27

-----  
Standard Lot#: 25414100A4  
Last Changed By: T. GALLEGOS #28847  
-----

OAS: W. LOOMIS #53545  
TUCSON PD  
-----

Operator: W. LOOMIS #53545  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:18:56
0.100 Cal Check	0.097	09:19:17
Air Blank	0.000	09:19:47

Preventative Maintenance  
Performed -- Timer Reset

|