



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TS4
2-9-2016


I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002760 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 2-5-16 Cal 1 0.096 AC Cal 2 0.095 AC

Between the dates of 1-8-16 and 2-5-16 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on .
The instrument was returned to service on .

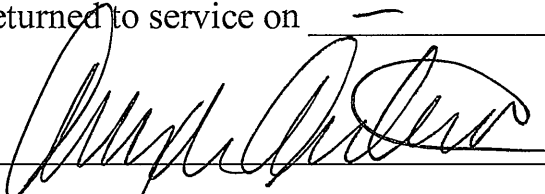

Amy L. Autret

2-9-16

ADPS QAS Permit # 33718

SUBSCRIBED AND SWORN to before me this 9th, day of February, 2016.



Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

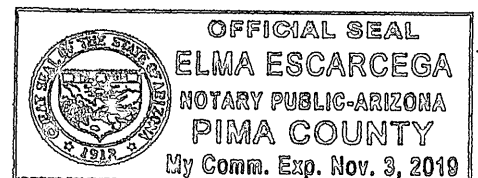


EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7
JH
TS4
2.9.2016

QA SPECIALIST Amy Autret 53546 AGENCY Tucson Police Department

DATE 2-5-16 TIME 0919

INTOXILYZER SERIAL # 80-002760 LOCATION ODE-Intox room

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
(✓) 2. Date check.


OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 095 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS

SIGNATURE

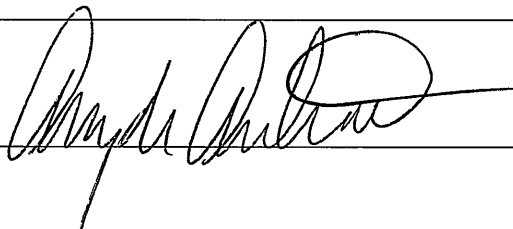

INTOXILYZER 8000

Location: TUCSON PD

Serial Number: 80-002760

Core Version: 8105.48

02/05/2016

09:20:45

Standard Lot#: 25414100A4

Last Changed By: T. GALLEGOS

#28847

QAS: A. AUTRET

#53546

TUCSON PD

Operator: A. AUTRET

#53546

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1987

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:21:48
Diagnostic Test	Pass	09:22:23
Air Blank	0.000	09:22:50
0.100 Cal Check	0.096	09:23:11
Air Blank	0.000	09:23:41
Subject Test	0.000	09:24:07
Air Blank	0.000	09:24:34
Wait	ABT*	
Air Blank	0.000	09:25:08

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2 of 7
JH
TS4

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002760
Core Version: 8105.48
02/05/2016

09:25:12

Standard Lot#: 25414100A4
Last Changed By: T. GALLEGOS #28847

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:26:06
Diagnostic Test	Pass	09:26:41
Air Blank	0.000	09:27:09
0.100 Cal Check	RFI*	09:27:24
Air Blank	IPS**	09:27:33

*RFI Detect

**Improper Sample

3 of 7
[Signature]
TS4

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002760
Core Version: 8105.48
02/05/2016

09:27:37

4 of 7
JLB
154

Standard Lot#: 25414100A4
Last Changed By: T. GALLEGOS #28847

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:28:37
Diagnostic Test	Pass	09:29:12
Air Blank	0.000	09:29:39
0.100 Cal Check	0.097	09:30:00
Air Blank	0.000	09:30:29
Subject Test	0.114*	09:33:37
Air Blank	0.000	09:34:10
Air Blank	0.000	09:34:37
Subject Test	XXX**	09:35:02
Air Blank	0.000	09:35:36

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON PD

Serial Number: 80-002760

Core Version: 8105.48

02/05/2016

09:35:54

Standard Lot#: 25414100A4

Last Changed By: T. GALLEGOS #28847

QAS: A. AUTRET #53546

TUCSON PD

Operator: A. AUTRET #53546

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:36:23
0.100 Cal Check	0.095	09:36:43
Air Blank	0.000	09:37:13

Preventative Maintenance

Performed -- Timer Reset

5 of 7
JLB
TBH

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY


STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Amy Autret 53546 AGENCY Tucson Police Department
DATE 2-5-16 TIME 0919
Intoxilyzer Serial # 80-002760 Location ODE - Intox room

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 096 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

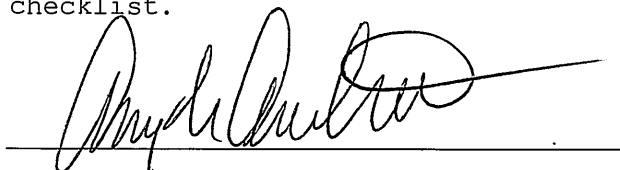

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002760
Core Version: 8105.48
02/05/2016 09:19:13

Standard Lot#: 25414100A4
Last Changed By: T. GALLEGOS #28847

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:19:41
0.100 Cal Check	0.096	09:20:02
Air Blank	0.000	09:20:31

Preventative Maintenance
Performed -- Timer Reset

7 of 7
[Signature]
Bh