



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002760 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TSy
4.7-2016

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3/31/2016 Cal 1 0.097 AC Cal 2 0.095 AC

Between the dates of 3/4/2016 and 3/31/2016 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on
The instrument was returned to service on

William L Loomis

William L. Loomis

4/4/2016


ADPS QAS Permit # 33720

SUBSCRIBED AND SWORN to before me this 4th day of April, 2016.

Elma Escarcega

Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

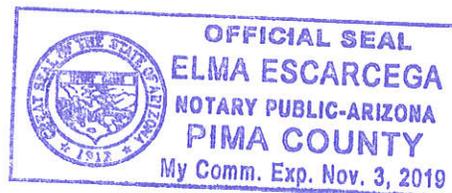

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7  
WW  
Apr 7. 2016

QA SPECIALIST W. LOOMIS 53545 AGENCY Tucson Police Department

DATE 3/31/2016 TIME 10:46

INTOXILYZER SERIAL # 80-002760 LOCATION ODE - SGT'S OFFICE

(  ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(  ) 1. Clock time check.

(  ) 2. Date check.

OPERATIONAL TESTS

(  ) 1. Alcohol-free subject test result 0. 000 AC.

(  ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed

(  ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.

(  ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 095 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE William Loomis

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002760  
Core Version: 8105.48  
03/31/2016 10:48:00

2 of 7  
WLL  
TS4

Standard Lot#: 25414100A4  
Last Changed By: T. GALLEGOS #28847

OAS: W. LOOMIS #53545  
TUCSON PD

Operator: W. LOOMIS #53545  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1987  
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:48:47
Diagnostic Test	Pass	10:49:24
Air Blank	0.000	10:49:51
0.100 Cal Check	0.097	10:50:12
Air Blank	0.000	10:50:41
Subject Test	0.000	10:50:58
Air Blank	0.000	10:51:26
Wait	ABT*	
Air Blank	0.000	10:52:00

\*Sequence Aborted

Not a Successfully  
Completed  
Test Sequence

3 of 7  
WU  
TSh

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002760  
Core Version: 8105.48  
03/31/2016 10:52:07

Standard Lot#: 25414100A4  
Last Changed By: T. GALLEGOS #28847

OAS: W. LOOMIS #53545  
TUCSON PD

Operator: W. LOOMIS #53545  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1987  
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:53:51
Diagnostic Test	Pass	10:54:28
Air Blank	0.000	10:54:56
0.100 Cal Check	RFI*	10:55:09
Air Blank	IPS**	10:55:17

\*RFI Detect  
\*\*Improper Sample

Not a Successfully  
Completed  
Test Sequence

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002760  
Core Version: 8105.48  
03/31/2016 10:55:32

4 of 7  
WLL  
PSG

Standard Lot#: 25414100A4  
Last Changed By: T. GALLEGOS #28847

OAS: W. LOOMIS #53545  
TUCSON PD

Operator: W. LOOMIS #53545  
TUCSON PD

Subject: SOAP, SOAP  
DOB: 11/11/1987  
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:56:27
Diagnostic Test	Pass	10:57:05
Air Blank	0.000	10:57:32
0.100 Cal Check	0.096	10:57:53
Air Blank	0.000	10:58:22
Subject Test	0.143*	11:01:30
Air Blank	0.000	11:02:03
Air Blank	0.000	11:02:31
Subject Test	XXX**	11:02:47
Air Blank	0.000	11:03:22

\*Deficient Sample  
\*\*Invalid Sample  
Begin new deprivation period

Not a Successfully  
Completed  
Test Sequence

5 of 7  
wll  
TS4

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002760  
Core Version: 8105.48  
03/31/2016 11:03:41

-----  
Standard Lot#: 25414100A4  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: W. LOOMIS #53545  
TUCSON PD  
-----

Operator: W. LOOMIS #53545  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	11:04:10
0.100 Cal Check	0.095	11:04:31
Air Blank	0.000	11:05:01

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD CALIBRATION CHECK PROCEDURE

6 of 7  
WLL  
TS4  
4.7.2016

QA Specialist W. LOOMIS 53545 AGENCY Tucson Police Department  
DATE 3/31/2016 TIME 10:46  
Intoxilyzer Serial # 80-002760 Location ODE-SGT'S OFFICE

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 097 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Willie Loomis

DPS Form Exh G-2 (Rev. 05-01)

7 of 7  
TSJ WLL

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002760  
Core Version: 8105.48  
03/31/2016  
10:46:18

Standard Lot#: 25414100A4  
Last Changed By: T. GALLEGOS #28847

OAS: W. LOOMIS #53545  
TUCSON PD

Operator: W. LOOMIS #53545  
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:46:47
0.100 Cal Check	0.097	10:47:08
Air Blank	0.000	10:47:38

Preventative Maintenance  
Performed - Timer Reset