



CERTIFICATION by CUSTODIAN of RECORDS  
Tucson Police Department Crime Laboratory

Ch  
2.20.08

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002803 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 2/14/08 Cal 1 0.101 AC Cal 2 0.100 AC

Between the dates of 1/16/08 and 2/14/08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
The instrument was returned to service on -.


Sarah Cione

Sarah L. Cione

2/15/08

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 15, day of February, 2008.

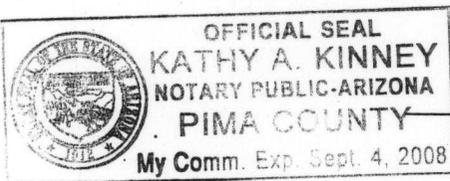

Kathy A Kinney  
Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1 of 9sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

Ch

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

2.20.08

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 2/14/08 TIME 9:18

INTOXILYZER SERIAL # 80-002803 LOCATION Team 2

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 100 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002803  
02/14/2008

09:19:37

2 of 9 sc

Ch

2-20-08

Standard Lot#: 732503I  
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:20:27
Diagnostic Test	Pass	09:20:55
Air Blank	0.000	09:21:22
0.100 Cal Check	0.100	09:21:41
Air Blank	0.000	09:22:11
Subject Test	0.000	09:22:27
Air Blank	0.000	09:22:54
Five Minute Wait	ABT*	
Air Blank	0.000	09:23:30

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002803  
02/14/2008

09:23:38

3 of 95c

OK

2-20-08

---

Standard Lot#: 732503I  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP,TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	09:24:30
Diagnostic Test	Pass	09:24:57
Air Blank	RFI*	09:25:06
Air Blank	IPS**	09:25:15

\*RFI Detect

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002803  
02/14/2008

09:25:21

4 of 95c  
Ln  
2-22-08

Standard Lot#: 732503I  
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP,TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:26:12
Diagnostic Test	Pass	09:26:40
Air Blank	0.000	09:27:07
0.100 Cal Check	0.101	09:27:26
Air Blank	0.000	09:27:55
Subject Test	0.162*	09:31:04
Air Blank	0.000	09:31:36
Air Blank	0.000	09:32:04
Subject Test	XXX**	09:32:27
Air Blank	0.000	09:33:01

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002803  
02/14/2008 09:33:17

5 of 9 sc

OK  
2-20-08

-----  
Standard Lot#: 732503I  
Last Changed By: S. CIONE #50616

-----  
QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:33:45
Cal Check	0.100	09:34:03
Air Blank	0.000	09:34:33

*Sarah Cione*

Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

6 of 9sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

Ch

2-20-08

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 2/14/08 TIME 9:18  
Intoxilyzer Serial # 80-002803 Location Team 2

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach  $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)


INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002803  
02/14/2008 09:18:11

7 of 9 SC  
CU  
2-22-08

-----  
Standard Lot#: 732503I  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:18:38
Cal Check	0.101	09:18:57
Air Blank	0.000	09:19:27

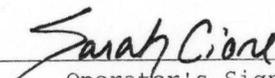  
-----  
Operator's Signature



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

8 of 95c

INTOXILYZER SERVICE RECORD  
FOR  
CALIBRATION GAS TANK CHANGE

lu  
2-20-08

Date: 2/14/08

QAS: Sarah Cione

Intoxilyzer Serial Number: 80-002803

New Lot Number: 732503 I

New Expiration Date: 11/26/2009

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."  
Select "D" for dry gas and enter the following information:  
  - Enter Target Value: 0.100
  - Enter Control Lot #: From the tank label
  - Enter Exp Date: From the tank label (mm/dd/yyyy)
  - Edit info Y/N? Y or N
  - Update Std Info? Y ENTER
  - "Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."  
Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 101 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove and sign printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002803  
02/14/2008 09:16:28

9 of 95c

Ch

2-22-08

-----  
Standard Lot#: 732503I  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:16:55
Cal Check	0.101	09:17:14
Air Blank	0.000	09:17:44

*Sarah Cione*

Operator's Signature